

# TEST DESIGN SPECIFICATION

KAU Hospital System (KHS)
Independence Verification and Validation

Version: 1.0.0

Date : 06/24/2024

Project Title: KHS Independent Verification and Validation		
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	



#### **Document Control**

Document Name KHS Test Design Specification	
Reference Number	KHS_TDS
Version	1.0.0
Project Code	SAADA_KHS
Status	In-use
Date Released	29 <sup>th</sup> of April, 2024

Name	Position	Signature
Prepared By: Danisha Azra	Test Lead	
Prepared By: Shahrul Amin	Test Automation Engineer	
Reviewed By: Nur Aisyah	Tester	
Approved By: Adiba Alya	Test Manager	
Approved By: Amylia Natasya	Test Analyst	

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## **Version History**

Version	Release Date	Section	Amendments
1.0.0	29/04/2024	All	Original Document

#### **Distribution List**

Version	Release Date	Controlled Copy No	Recipient Name	Department	Issue Date	Return Date
1.0.0	29/04/2024	01		Testerontop Test Team	29/04/2024	

Date:06/04/2024

Test Plan ID: KHS\_TDS\_1\_1.0.0



#### **Table of Content**

Document Control	
Version History	3
Distribution List	3
Table of Content	4
1.0 Introduction	6
1.1 Purpose	6
1.2 Scope	6
1.3 References	6
2.0 Test Design	7
2.1 Features to be Tested	7
2.2 Approach Refinements	7
2.2.1 F001 Add Doctor	8
2.2.1.1 Equivalence Partitioning	8
2.2.1.2 Boundary Value Analysis	11
2.2.1.3 Decision Table Testing	12
2.2.1.4 Use Case Testing	13
2.2.2 F002 Add Patient	15
2.2.2.1 Equivalence Partitioning	15
2.2.2.2 Boundary Value Analysis	19
2.2.2.3 Decision Table Testing	21
2.2.2.4 Use Case Testing	22
2.2.3 F003 Add Medicine	24
2.2.3.1 Equivalence Partitioning	24
2.2.3.2 Boundary Value Analysis	28
2.2.3.3 Decision Table Testing	30
2.2.3.4 Use Case Testing	31
2.2.4 F004 Add Nurse	33
2.2.4.1 Equivalence Partitioning	33
2.2.4.2 Boundary Value Analysis	36
2.2.4.3 Decision Table Testing	37
2.2.4.4 Use Case Testing	38
2.2.5 F005 Add Room	40
2.2.5.1 Equivalence Partitioning	40
2.2.5.2 Boundary Value Analysis	44
2.2.5.3 Decision Table Testing	47

Date:06/04/2024

Test Plan ID: KHS\_TDS\_1\_1.0.0



2.2.5.4 Use Case Testing	48
2.3 Test Identification	53
2.3.1 Add Doctor	53
2.3.2 Add Patient	6 <sup>2</sup>
2.3.3 Add Medicine	65
2.3.4 Add Nurse	75
2.3.5 Add Room	
2.4 Features Pass/Fail Criteria	93
2.5 Test Deliverables	QZ

Project Title: KHS Independen	A TERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE 250

#### 1.0 Introduction

#### 1.1 Purpose

This test design specification support the following objectives:

- 1. To identify the test items covered
- 2. To further refine the test approach from Test Plan
- 3. To define the sources of the information used to prepare the plan
- 4. To identify the general cases of test to be executed

#### 1.2 Scope

This test design specification covers Session Verification features developed in KHS.

#### 1.3 References

The following IEEE standards have been referenced in preparation of this document:

- 1. IEEE 829-1998 Standard for Software and System Test Documentation
- 2. IEEE 829-2008 Standards for Software and System Test Documentation

The following documents provide the test basis for this test design:

- 1. KHS Software Requirement Specifications (KHS SRS 1.0)
- 2. KHS System Design Specifications (KHS\_SDS\_1.0)

The following GitHub repositories provide the code basis for this test design:

 rahmaAlghalib. (2022, Aug 8). GitHub - rahmaAlghalib/KAU-Hospital-System: Java Hospital management system. GitHub.

https://github.com/rahmaAlghalib/KAU-Hospital-System

Project Title: KHS Independent	TERON A	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	A TO S TO

#### 2.0 Test Design

#### 2.1 Features to be Tested

The following table contains the features to be tested based on KHS Software Requirement Specifications (SRS) for this iteration. Listed together are the function ID and its corresponding functions and estimated risk level.

Function Id	Functions	Risk Level
F001	Add Doctor	High
F002	Add Patient	High
F003	Add Medicine	Medium
F004	Add Nurse	Medium
F005	Add Room	Medium

**Table 2.1 Features To Be Tested** 

#### 2.2 Approach Refinements

As discussed in Test Plan (KHS\_TP\_1.0.0), test cases for KHS will be designed using four techniques which are covered in this section. Techniques are applied based on suitability of techniques according to the nature of features.

The techniques that will be applied are as follows:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

Nature of feature as below:

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.1 F001 Add Doctor

Four techniques are identified for designing test for this feature, including:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

#### Equivalence partitioning is used to

#### 2.2.1.1 Equivalence Partitioning

### Equivalence classes for id

D1: 5000 ≤ doctorld ≤ 7999

D2: doctorld < 5000 D3: doctorld > 7999

#### **Equivalence classes for gender**

G1: gender  $\in$  {'F','M'} G2: gender  $\notin$  {'F','M'}

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



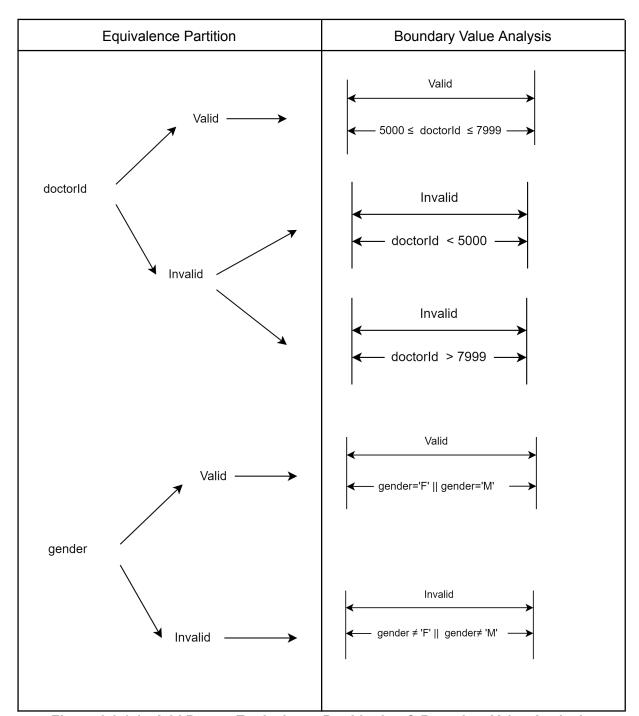


Figure 2.2.1.1a Add Doctor Equivalence Partitioning & Boundary Value Analysis

Date:06/04/2024 Test Plan ID:



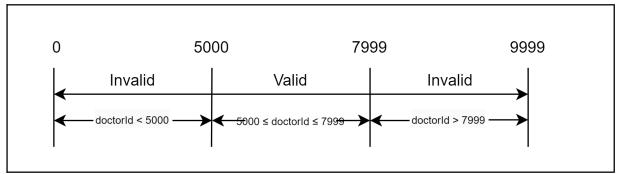


Figure 2.2.1.1b Add Doctor Equivalence Partitions


Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-01-001	doctorld< 5000	TCOV-01-001	doctorId < 5000 (Invalid)	999
TCON-01-002	5000 ≤ doctorId ≤ 7999	TCOV-01-002	5000 ≤ doctorId ≤ 7999 (Valid)	1234
TCON-01-003	doctorId> 7999	TCOV-01-003	doctorId > 7999 (Invalid)	8888
TCON-01-004	Valid doctorId	TCOV-01-004	Valid doctorId (valid output)	1234
TCON-01-005	Invalid doctorld	TCOV-01-005	Invalid doctorld (Invalid output)	12346
TCON-01-006	gender ∉ {'F','M'}	TCOV-01-006	gender ∉ {'F','M'} (invalid)	В
TCON-01-007	gender ∈ {'F','M'}	TCOV-01-007	gender ∈ {'F','M'} (valid)	F
TCON-01-008	Valid gender	TCOV-01-008	valid gender (valid output)	F
TCON-01-009	Invalid gender	TCOV-01-009	Invalid gender (invalid output)	В

Table 2.2.1.1 Add Doctor Equivalence Partitioning Test Conditions & Coverage

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



## 2.2.1.2 Boundary Value Analysis

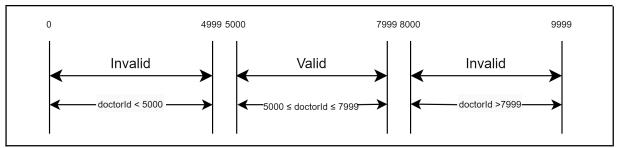

2.2.1.2 Add Doctor Boundary Value Analysis

EP1: 5000 ≤ doctorId ≤ 7999 EP2: 0 ≤ doctorId ≤ 4999 EP3: 8000 ≤ doctorId ≤ 9999

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-01-010	doctorId = 4999	TCOV-01-010	doctorld = 4999 (Invalid)	4999
TCON-01-011	doctorId = 5000	TCOV-01-011	doctorld = 5000 (Valid)	5000
TCON-01-012	doctorld = 5001	TCOV-01-012	doctorld = 5001 (Valid)	5001
TCON-01-013	doctorId = 7998	TCOV-01-013	doctorld = 7998 (Valid)	7998
TCON-01-014	doctorId = 7999	TCOV-01-014	doctorld = 7999 (Valid)	7999
TCON-01-015	doctorId = 8001	TCOV-01-015	doctorId = 8001(Invalid)	8001
TCON-01-004	Valid doctorId	TCOV-01-004	Valid doctorId (Valid output)	5000
TCON-01-005	Invalid doctorld	TCOV-01-005	Invalid doctorld (Invalid output)	8001

Table 2.2.1.2 Add Doctor Boundary Value Analysis Test Conditions & Coverage

Project Title: KHS Independent Ver	rification and Validation
------------------------------------	---------------------------

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.2.1.3 Decision Table Testing

Test Condition ID	Test Condition	
TCON-01-016	doctorId	
TCON-01-017	gender = 'F'	
TCON-01-018	gender = 'M'	
TCON-01-019	Display "Invalid doctor ID"	
TCON-01-020	Display "Invalid gender"	
TCON-01-021	Display "Doctor successfully added"	

Table 2.2.1.3a Add Doctor Decision Table Test Conditions

Condition	1	2	3	4	5
doctorId	5000 ≤ doctorId ≤ 7999	5000 ≤ doctorId ≤ 7999	5000 ≤ doctorId ≤ 7999	doctorId < 5000	doctorId > 7999
gender = 'F'	Υ	N	N	-	-
gender = 'M'	N	Υ	N	-	-
Output					
Display "Invalid ID"				X	Х
Display "Invalid gender"			X		
Display "Doctor successfully added"	x	X			

**Table 2.2.1.3b Add Doctor Decision Table** 

Project Title: KHS Independent	4 ERONA O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE 202

## 2.2.1.4 Use Case Testing

Use Case ID	UC001			
Use Case	F001 Add Doctor			
Purpose	To add o	To add doctor for database		
Actor	Hospital	Administrator		
Description	The hospital administrator starts by providing the information needed to add nurses into the database. If all valid requirements are input, the doctor will be added successfully into the database.			
Precondition	System	prompt options		
Post-condition	Doctor s	successfully added to the database		
Scenario Name	Step Action			
Main Flow	1	System choose option '1'		
	2	User input doctorld		
	3	User input gender		
	4	System verify valid doctorld		
	5	System verify gender		
	6	System add doctor to the database		
	7	System displays message "Doctor successfully added"		
Alternate Flow -	2.1	User input invalid doctorId		
Invalid doctorld	2.2	System display error message "Invalid doctor ID"		
	2.3	The use case end		
Alternate Flow -	3.1	User input invalid doctorId		
Invalid Gender	3.2	System displays error message "Invalid gender"		
	3.3	The use case end		

Project Title: KHS Independent	4 ERONA	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE SOL

Rules	None
Note	None

Table 2.2.1.4a Add Doctor Use Case

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Test Data
TCON-01-022	Main flow	TCOV-01-022	Main Flow	Valid doctorId = 5555 Valid Gender = 'F'
TCON-01-023	Alternate Flow - Invalid doctorId	TCOV-01-023	Alternate Flow - Invalid doctorId	Invalid doctorId = 1234
TCON-01-024	Alternate Flow - Invalid gender	TCOV-01-024	Alternate Flow - Invalid gender	Valid doctorld = 5555 Invalid Gender = 'B'

Table 2.2.1.4b Add Doctor Use Case Test Conditions & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.2 F002 Add Patient

Four techniques are identified for designing test for this feature, including:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

#### Equivalence partitioning is used to

#### 2.2.2.1 Equivalence Partitioning

#### **Equivalence classes for patientID**

A1 : 1000 ≤ patientID ≤ 3999

A2 : patientID < 1000 A3 : patientID > 3999

#### **Equivalence classes for illness**

C1 : illnessLength > 0 C2 : illnessLength  $\leq$  0

Date:06/04/2024 Test Plan ID:



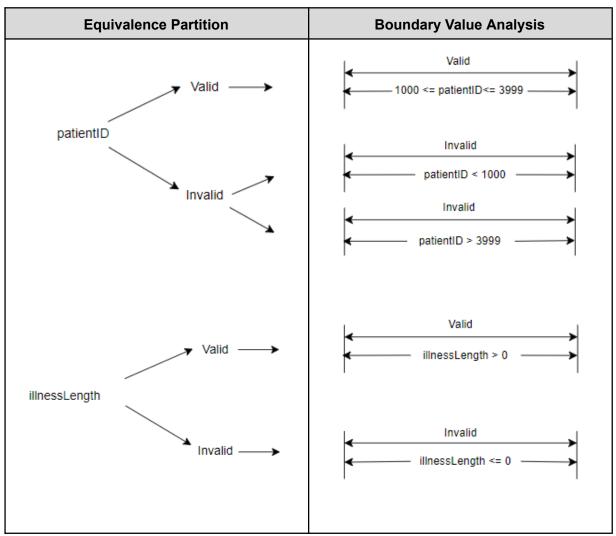


Figure 2.2.2.1a Add Doctor Equivalence Partitioning & Boundary Value Analysis

Project Title: KHS Independent Verification and Validation
--

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



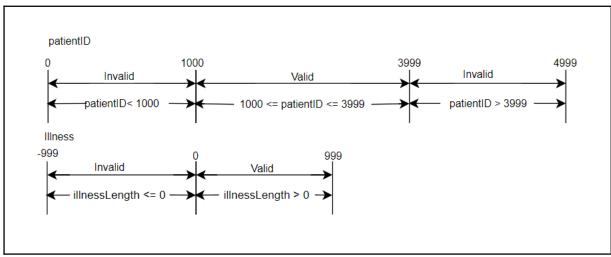


Figure 2.2.2.1b Add Patient Equivalence Partitioning

Date:06/04/2024 Test Plan ID:



Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-02-001	1000 ≤ patientID ≤ 3999	TCOV-02-001	1000 ≤ patientID ≤ 3999 (Valid)	3870
TCON-02-002	patientID < 1000	TCOV-02-002	patientID < 1000 (Invalid)	600
TCON-02-003	patientID > 3999	TCOV-02-003	patientID > 3999 (Invalid)	5678
TCON-02-004	Valid patientID	TCOV-02-004	Valid patientID (valid output)	1234
TCON-02-005	Invalid patientID	TCOV-02-005	Invalid patientID (invalid output)	123
TCON-02-006	illnessLength > 0	TCOV-02-006	illnessLength > 0 (Valid)	'Flu'
TCON-02-007	illnessLength <= 0	TCOV-02-007	illnessLength <= 0 (Invalid)	σ
TCON-02-008	Valid illnessLength	TCOV-02-008	Valid illnessLength (valid output)	'Fever'
TCON-02-009	Invalid illnessLength	TCOV-02-009	Invalid illnessLength (Invalid output)	ø

Table 2.2.2.1 Add Patient Equivalence Partitioning Test Condition & Coverage

Date:06/04/2024 Test Plan ID:

Test Plan ID:
KHS\_TDS\_1\_1.0.0

## 2.2.2.2 Boundary Value Analysis

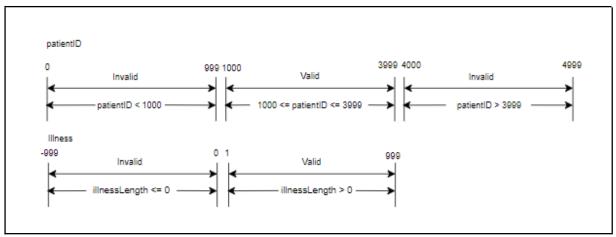


Figure 2.2.2.2 Add Patient Boundary Value Analysis

EP1 : 1000 ≤ patientID≤ 3999

EP2 :patientID < 1000 EP3 :patientID > 3999

EP4 : illnessLength > 0 EP5 : illnessLength ≤ 0

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-02-010	patientID = 999	TCOV-02-010	patientID = 999 (Invalid)	999
TCON-02-011	patientID = 1000	TCOV-02-011	patientID= 1000 (valid)	1000
TCON-02-012	patientID = 1001	TCOV-02-012	patientID = 1001 (valid)	1001
TCON-02-013	patientID = 3998	TCOV-02-013	patientID= 3998 (valid)	3998

Date:06/04/2024 Test Plan ID:



TCON-02-014	patientID= 3999	TCOV-02-014	patientID = 3999 (valid)	3999
TCON-02-015	N-02-015 patientID= 4000 TCOV-02-015		patientID = 4000 (invalid)	4000
TCON-02-004	Valid patientID	TCOV-02-016	Valid output	1234
TCON-02-005	Invalid patientID	TCOV-02-017	Invalid output	4000
TCON-02-018	illnessLength = 0	TCOV-02-018	illnessLength = 0 (invalid)	σ
TCON-02-019	illnessLength = 1	TCOV-02-019	illnessLength =1 (valid)	'F'
TCON-02-008	Valid illnessLength	TCOV-02-020	Valid output	'Cold'
TCON-02-009	Invalid illnessLength	TCOV-02-021	Invalid output	tı

Project Title: KHS	S Independent	Verification and	d Validation
T TO COL TILIC. IN IN	, illiacpellacili	v Ci ilicationi ani	a vanuation

Date:06/04/2024 Test Plan I

Test Plan ID: KHS\_TDS\_1\_1.0.0



# 2.2.2.3 Decision Table Testing

Test Condition ID	Test Condition
TCON-02-022	patientID
TCON-02-023	IllnessLength
TCON-02-024	Displays an error message "Invalid patient ID"
TCON-02-025	Displays an error message "Illness cannot be null"
TCON-02-026	Displays message "Patient successfully added"

Condition	1	2	3	4
patientID	<1000	1000 ≤ ID ≤ 3999	1000 ≤ ID ≤ 3999	ID > 3999
Illness Length	-	≤ 0	> 0	-
Output				
Displays an error message "Error: Invalid patientID"	×			х
Displays an error message "Error: Illness cannot be empty."		×		
Displays message "Patient successfully added"			Х	

Project Title: KHS Independent	SERONA.	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SUDCE SOL



## 2.2.2.4 Use Case Testing

Use Case ID	UC002		
Use Case	F002	Add Patient	
Purpose	To add	d a new patient to the hospital database	
Actor	Hospi	tal Administrator	
Description	The hospital administrator provides the information needed to add a new patient to the hospital database. The system verifies the provided information and, if valid, adds the new patient to the database and updates the system accordingly.		
Precondition	Syste	m prompt options	
Post-condition	Patier	nt successfully added to the database	
Scenario Name	Step	Action	
Main Flow	1	User choose option '4'	
	2	User input patientID	
	3	User input Illness	
	4	System verify valid patientID	
	5	System verify valid Illness	
	6	System adds new patient to the database	
	7	System displays message "Patient successfully added"	
Alternate Flow -	4.1	User input invalid patientID	
Invalid ID	4.2	System displays error message"Invalid patient ID"	
	4.3	The use case end	
Alternate Flow -	5.1	User input invalid Illness	
Invalid illness	5.2	System displays an error message "Illness cannot be empty."	
	5.3	The use case end	

Project Title: KHS Independent	STERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	Since zor

Rules	None
Note	None

**Table 2.2.2.4a Session Verification Use Case** 

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Test Data
TCON-02-02 7	Main Flow	TCOV-02-027	Main Flow	Valid patientID = 1234 Valid Illness = 'FLu'
TCON-02-02 8	Alternate Flow - Invalid patientID	TCOV-02-028	Alternate Flow - Invalid patientID	Invalid patientID = 124 Valid Illness = 'Flu'
TCON-02-02 9	Alternate Flow - Invalid Illness	TCOV-029	Alternate Flow - Invalid Illness	Valid patientID = 1234 Invalid Illness = "

Table 2.2.2.4b Session Verification Use Case Test Conditions & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.3 F003 Add Medicine

Four techniques are identified for designing test for this feature, including:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

Equivalence partitioning is used to split medicineCode and price into possible partitions for testing. Boundary value analysis will test the boundary values of the medicineCode and price input. Decision table is used to cover major branches of decision

#### 2.2.3.1 Equivalence Partitioning

#### **Equivalence classes for medicine code**

A1:  $4000 \le \text{medicineCode} \le 4514$ 

A2: medicineCode < 4000 A3: medicine Code > 4514

#### **Equivalence classes for price**

B1:  $1.00 \le price \le 50.00$ 

B2: price < 1.00 B3: price > 50.00

Date:06/04/2024 | Test Plan ID:



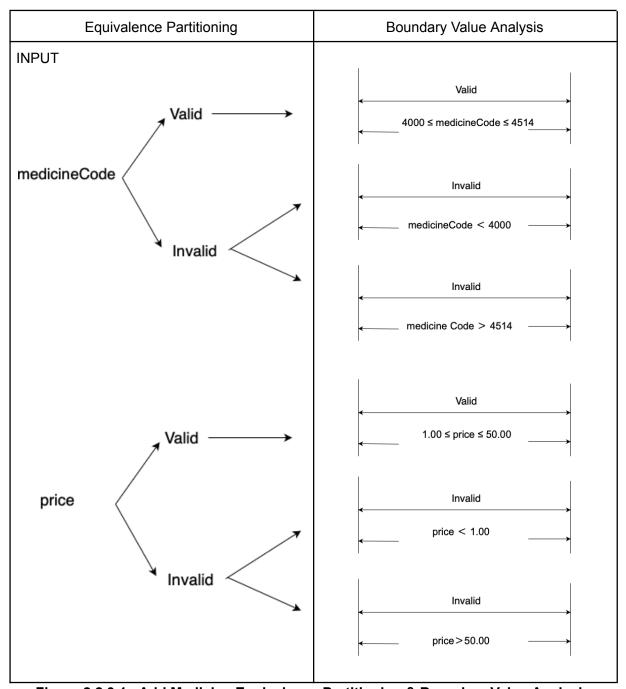


Figure 2.2.3.1a Add Medicine Equivalence Partitioning & Boundary Value Analysis

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



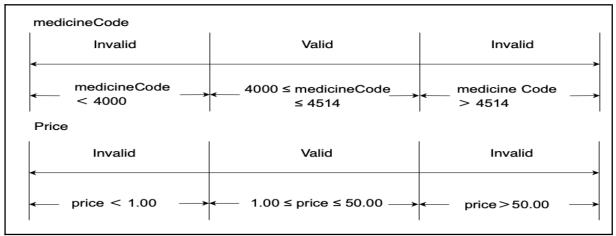


Figure 2.2.3.1b Add Medicine eEquivalence Partitioning

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-03-001	4000 ≤ medicineCode ≤ 4514	TCOV-03-001	4000 ≤ medicineCode ≤ 4514 (Valid)	4012
TCON-03-002	medicineCode < 4000			3890
TCON-03-003	medicine Code TCOV-03-003 > 4514		medicine Code > 4514 (Invalid)	4900
TCON-03-004	1.00 ≤ price ≤ 50.00			20.00
TCON-03-005	price < 1.00 TCOV-03-005		price < 1.00 (Invalid)	0.00
TCON-03-006	price>50.00	TCOV-03-006	Price>50.00 (Invalid)	1000.00

Project Title: KHS Independent	Verification and Validation	A ERON O
Date:06/04/2024	Test Plan ID: KHS TDS 1 1.0.0	SINCE 202

TCON-03-007	Valid medicineCode	TCOV-03-007	Valid Output	4210
TCON-03-008	Invalid medicineCode	TCOV-03-008	Invalid Output	4999
TCON-03-009	Valid price	TCOV-03-009	Valid Output	30.00
TCON-03-010	Invalid price	TCOV-03-010	Invalid Output	1000.00


Table 2.2.3.1 Add Medicine Equivalence Partitioning Test Condition & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.3.2 Boundary Value Analysis

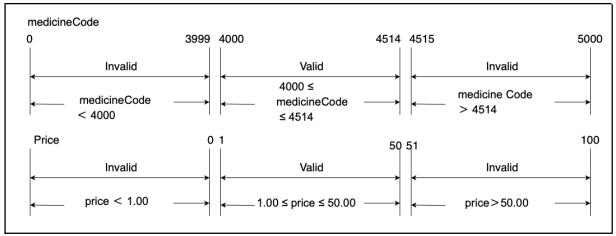

2.2.3.2 Input Boundary Value Analysis

EP1: 4000 ≤ medicineCode ≤ 4514

EP2: medicineCode < 4000 EP3: medicineCode > 4514

EP4:  $1.00 \le price \le 50.00$ 

EP5: price < 1.00 EP6: price > 50.00

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-03-011	medicineCode = 3999	TCOV-03-011	medicineCode = 3999 (Invalid)	3999
TCON-03-012	medicineCode = 4000	TCOV-03-012	medicineCode = 4000 (Valid)	4000
TCON-03-013	medicine Code = 4001	TCOV-03-013	medicine Code = 4001 (Valid)	4001
TCON-03-014	medicineCode = 4513	TCOV-03-014	medicineCode = 4513 (Valid)	4513

Date:06/04/2024 Test Plan ID:



TCON-03-015	medicineCode = 4514	TCOV-03-015	medicineCode = 4514 (Valid)	4514
TCON-03-016	medicineCode = 4515	TCOV-03-016	medicineCode = 4515 (Invalid)	4515
TCON-03-017	price = 0.00	TCOV-03-017	COV-03-017 price = 0.00 (Invalid)	
TCON-03-018	price = 1.00	TCOV-03-018	price = 1.00 (Valid)	1.00
TCON-03-019	price = 2.00	TCOV-03-019	price = 2.00 (Valid)	2.00
TCON-03-020	price =49.00	TCOV-03-020	price = 49.00 (Valid)	49.00
TCON-03-021	price = 50.00	TCOV-03-021	price = 50.00 (Valid)	50.00
TCON-03-022	price = 51.00	TCOV-03-022	price = 51.00 (Invalid)	51.00
TCON-03-007	Valid medicineCode	TCOV-03-007	Valid Output	4210
TCON-03-008	Invalid medicineCode	TCOV-03-008	Invalid Output	4999
TCON-03-009	Valid price	TCOV-03-009	Valid Output	30.00
TCON-03-010	Invalid price	TCOV-03-010	Invalid Output	1000.00

Table 2.2.3.2 Input Boundary Value Analysis Test Condition & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.2.3.3 Decision Table Testing

Test Condition ID	Test Condition	
TCON-03-023	medicineCode (mc)	
TCON-03-024	Price (p)	
TCON-03-025	Display "invalid medicine code"	
TCON-03-026	Display "invalid price"	
TCON-03-027	Add medicine into database	
TCON-03-028	Display "Medicine successfully added"	

Table 2.2.3.3a Add Medicine Decision Table Test Condition

Condition	TCOV-03-023	TCOV-03-024	TCOV-03-025	TCOV-03-026	TCOV-03-027
Medicine Code (mc)	mc<4000	4000≤ mc≤4514	4000≤ mc≤4514	4000≤ mc≤4514	mc>4514
Price(p)	-	p < 1.00	1.00 ≤p≤ 50.00	p > 50.00	-
Output					
Display "invalid medicine code"	*				×
Display "invalid price"		×		×	
Add medicine into database			*		
Display "Medicine successfu			*		

Project Title: KHS Independent	STERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE ZON

Ily added"
------------

**Table 2.2.3.3b Add Medicine Decision Table** 

# 2.2.3.4 Use Case Testing

Use Case ID	UC00	3		
Use Case	F003 Add Medicine			
Purpose	To add	To add new medicine into database		
Actor	Hospi	tal Administrator		
Description	The hospital administrator starts by providing the information needed to add new medicine into the database. If all valid requirements are input, the nurse will be added successfully into the database.			
Precondition	System prompt options			
Post-condition	Medicine successfully added into database			
Scenario Name	Step Action			
Main Flow	1 User choose option '5'			
	2	User enter medicine Code		
	3 User enter price			
	4 System validates medicine code			
	5	System validate price		
	6 System add nurse into database			
	7	System display "Medicine added successfully"		
Alternate Flow -	2.1	User Input invalidmedicineCode		
Invalid medicineCode	2.2	System display error messages "Invalid medicineCode"		
	2.3	The use case end		
Alternate Flow - Invalid price	4.1	4.1 User Input invalid price		

Project Title: KHS Independent	STERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE SOL

	4.2	System display error messages "Invalid price"
4.3 The use case end		The use case end
Rules	None	
Note	None	

Table 2.2.3.4a Add Medicine Use Case

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Test Data
TCON-03-02 9	Main Flow	TCOV-03-028	Main Flow	Name = paracetamol Valid medicine Code = 4102 Valid Price = 20.00
TCON-03-03 0	Alternate Flow - Invalid medicine Code	TCOV-03-029	Alternate Flow - Invalid medicine Code	Name = paracetamol Invalid medicine Code = 6234
TCON-03-03 1	Alternate Flow - Invalid Price	TCOV-03-030	Alternate Flow - Invalid Price	Name = paracetamol Valid medicine Code = 4001 Valid Price = 100.00

Table 2.2.3.4b Add Medicine Use Case Test Conditions & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.4 F004 Add Nurse

Four techniques are identified for designing test for this feature, including:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

Equivalence partitioning is used to split Id and gender into possible partitions for testing. Boundary value analysis will test the boundary values of the Id. Decision table is used to cover major branches of decision

#### 2.2.4.1 Equivalence Partitioning

#### **Equivalence classes for Id**

A1:  $8000 \le \text{nurseld} \le 9999$ 

A2: nurseld < 8000 A3: nurseld > 9999

#### Equivalence classes for gender

C1: gender ∈ {'M', 'F'} C2: gender ∉ {'M', 'F'}

Date:06/04/2024 Test Plan ID:



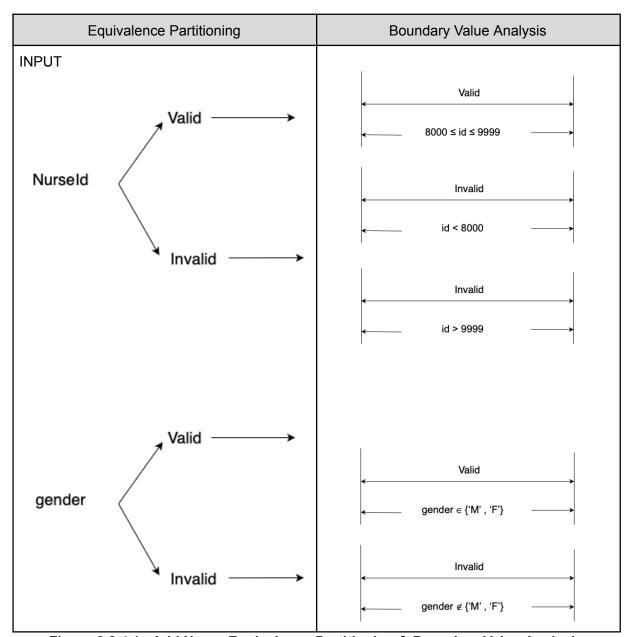


Figure 2.2.4.1a Add Nurse Equivalence Partitioning & Boundary Value Analysis

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



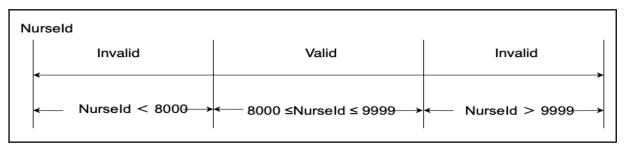


Figure 2.2.4.1b Add Nurse Equivalence Partitioning

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-04-001	8000 ≤ nurseId ≤ 9999	TCOV-04-001	8000 ≤ nurseld ≤9999 (Valid)	8125
TCON-04-002	nurseld < 8000	TCOV-04-002	nurseld < 8000 (Invalid)	6222
TCON-04-003	nurseld >9999	TCOV-04-003	nurseld >9999 (Invalid)	10000
TCON-04-004	gender ∈ {'M', 'F'}	TCOV-04-004	gender ∈ {'M', 'F'} (Valid)	М
TCON-04-005	gender ∉ {'M', 'F'}	TCOV-04-005	gender ∉ {'M','F'} (Invalid)	G
TCON-04-006	Valid nurseld	TCOV-04-006	Valid Output	8798
TCON-04-007	Invalid nurseld	TCOV-04-007	Invalid Output	4999
TCON-04-008	Valid gender	TCOV-04-008	Valid Output	F
TCON-04-009	Invalid gender	TCOV-04-009	Invalid Output	Α


Table 2.2.4.1 Add Medicine Equivalence Partitioning Test Condition & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.4.2 Boundary Value Analysis

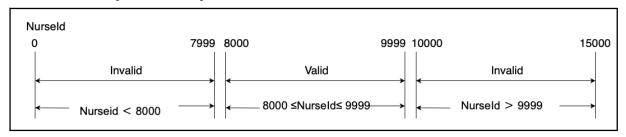

2.2.4.2 Input Boundary Value Analysis

EP1: 8000 ≤ nurseld ≤ 9000

EP2: nurseld < 8000 EP3: nurseld > 9999

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-04-010	nurseld = 7999	TCOV-04-010	nurseld = 7999 (Invalid)	7999
TCON-04-011	nurseld = 8000	TCOV-04-011	nurseld = 8000 (Valid)	8000
TCON-04-012	nurseld = 8001	TCOV-04-012	nurseld = 8001 (Valid)	8001
TCON-04-013	nurseld = 9998	TCOV-04-013	nurseld = 9998 (Valid)	9998
TCON-04-014	nurseld = 9999	TCOV-04-014	nurseld = 9999 (Valid)	9999
TCON-04-015	nurseld = 10000	TCOV-04-015	nurseld = 10000 (Invalid)	10000
TCON-04-008	Valid nurseld	TCOV-04-008	Valid Output	8798
TCON-04-009	Invalid nurseld	TCOV-04-009	Invalid Output	4999

Table 2.2.4.2 Input Boundary Value Analysis Test Condition & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.2.4.3 Decision Table Testing

Test Condition ID	Test Condition
TCON-04-016	nurseld
TCON-04-017	gender = 'M'
TCON-04-018	gender = 'F'
TCON-04-019	Display "Invalid nurse ID"
TCON-04-020	Display "Invalid gender"
TCON-04-021	Display "Nurse successfully added"

**Table 2.2.4.3a Add Nurse Decision Table Test Condition** 

Condition	1	2	3	5	7
nurseld	nurseld < 8000	8000 ≤ nurseId ≤ 9000	8000 ≤ nurseId ≤ 9000	8000 ≤ nurseId ≤ 9000	nurseld > 9999
gender = 'M'	-	Y	-	N	-
gender = 'F'	-	-	Υ	N	-
Output					
Print "Invalid nurse ID"	×				×
Print "Invalid gender"				×	
Display "Nurse successfully added"		×	×		

**Table 2.2.4.3b Add Nurse Decision Table** 

Project Title: KHS Independent	TERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE 252

# 2.2.4.4 Use Case Testing

Use Case ID	UC00	UC004		
Use Case	F004 Add Nurse			
Purpose	To add	To add nurse for database		
Actor	Hospi	tal Administrator		
Description	The hospital administrator starts by providing the information needed to add nurses into the database. If all valid requirements are input, the nurse will be added successfully into the database.			
Precondition	Hospi	tal administrators enter command 'Add_Nurse'		
Post-condition	Nurse	successfully added into database		
Scenario Name	Step	Action		
Main Flow	1	User choose option '3'		
	2	User input nurseld		
	3	User input gender		
	4 System verify valid nurseld			
	5 System verify gender			
	6 System add nurse into database			
	7	System display messages "Nurse successfully added"		
Alternate Flow -	2.1	User Input invalid nurseld		
Invalid Id	2.2	System display error messages "Invalid nurse ID"		
	2.3	The use case end		
Alternate Flow -	3.1	User Input invalid gender		
Invalid Gender	3.2	System display error messages "Invalid gender"		
	3.3	The use case end		

Project Title: KHS Independent	4 ERONA	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE 202

Rules	None
Note	None

Table 2.2.4.4a Add Nurse Use Case

Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Test Data
TCON-04-022	Main Flow	TCOV-04-022	Main Flow	Valid Id= 8203 Valid Gender = 'M'
TCON-04-023	Alternate Flow - Invalid Id	TCOV-04-023	Alternate Flow - Invalid Id	Invalid Id= 2309
TCON-04-024	Alternate Flow - Invalid Gender	TCOV-04-024	Alternate Flow - Invalid Gender	Valid Id= 8203 Invalid Gender = 'A'

Table 2.2.4.4b Add Nurse Use Case Test Conditions & Coverage

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



#### 2.2.5 F005 Add Room

Four techniques are identified for designing test for this feature, including:

- i. Equivalent partitioning
- ii. Boundary value analysis
- iii. Decision table testing
- Iv. Use Case Testing

### Equivalence partitioning is used to

### 2.2.5.1 Equivalence Partitioning

### Equivalence classes for roomNo

R1: 1 ≤ roomNo ≤ 399

R2: roomNo < 1 R3: roomNo > 399

#### Equivalence classes for charges


C1: 100 ≤ charges ≤ 2000

C2: charges < 100 C3: charges > 2000

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



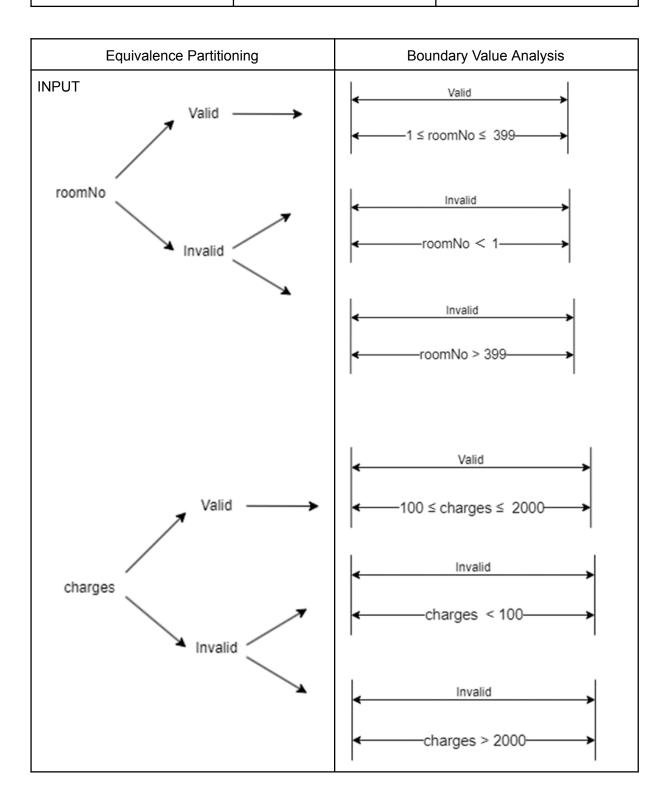

Project Title: KHS Independent	STERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE ZOZ

Figure 2.2.5.1a Add Room Equivalence Partitioning & Boundary Value Analysis

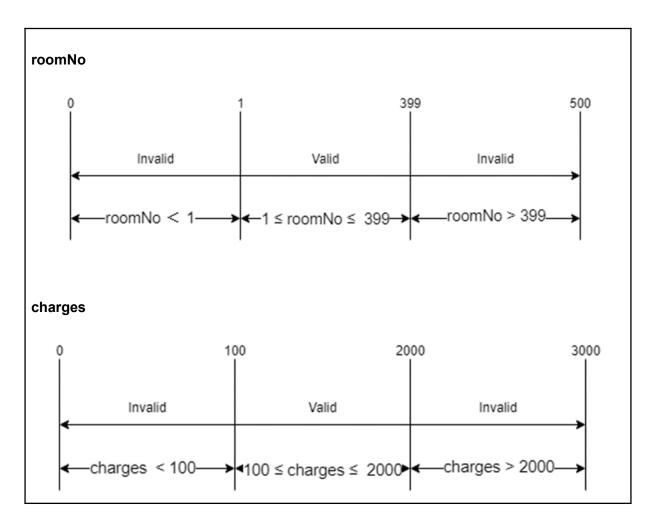


Figure 2.2.5.1b Add Room Equivalence Partitions

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0




Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-05-001	roomNo < 1	TCOV-05-001	roomNo < 1 (Invalid)	000
TCON-05-002	1 ≤ roomNo ≤ 399	TCOV-05-002	1 ≤ roomNo ≤ 399 (Valid)	101
TCON-05-003	roomNo > 399	TCOV-05-003	roomNo > 399 (Invalid)	443
TCON-05-004	Valid roomNo	TCOV-05-004	Valid roomNo (Valid output)	101
TCON-05-005	Invalid roomNo	TCOV-05-005	Invalid roomNo (Invalid output)	000
TCON-05-006	charges < 100	TCOV-05-006	charges < 100 (Invalid)	55
TCON-05-007	100 ≤ charges ≤ 2000	TCOV-05-007	100 ≤ charges ≤ 2000 (Valid)	900.0
TCON-05-008	charges > 2000	TCOV-05-008	charges > 2000 (invalid)	2500.00
TCON-05-009	Valid charges	TCOV-05-009	Valid charges (Valid output)	900.0
TCON-05-010	Invalid charges	TCOV-05-010	Invalid charges (Invalid output)	2500.0

Table 2.2.5.1c Add Room Equivalence Partitioning Test Conditions & Coverage

KHS\_TDS\_1\_1.0.0



## 2.2.5.2 Boundary Value Analysis

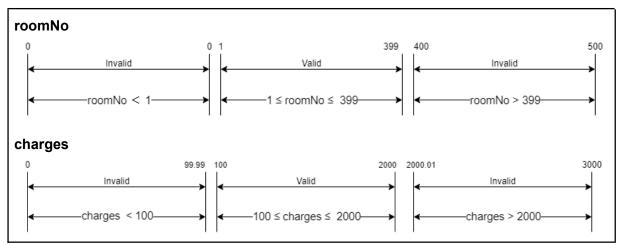

2.2.5.2a Input Boundary Value Analysis

EP1: 1 ≤ roomNo ≤ 399

EP2: roomNo < 1

EP3:  $400 \le \text{roomNo} \le 500$ 

EP4:  $100 \le \text{charges} \le 2000$ EP5:  $0 \le \text{charges} \le 99.99$ EP6:  $2000.01 \le \text{charges} \le 3000$ 

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Example of Test Data
TCON-05-011	roomNo = 000	TCOV-05-011	roomNo = 000 (Invalid)	000
TCON-05-012	roomNo = 001	TCOV-05-012	roomNo = 001 (Valid)	001
TCON-05-013	roomNo = 002	TCOV-05-013	roomNo = 002 (Valid)	002
TCON-05-014	roomNo = 398	TCOV-05-014	roomNo = 398	398
TCON-05-015	roomNo = 399	TCOV-05-015	roomNo = 399 (Valid)	399
TCON-05-016	roomNo = 400	TCOV-05-016	roomNo = 400 (Invalid)	400
TCON-05-017	charges = 99.99	TCOV-05-017	charges = 99.99 (Invalid)	99.99
TCON-05-018	charges = 100.00	TCOV-05-018	charges = 100 (Valid)	100
TCON-05-019	charges = 100.01	TCOV-05-019	charges = 100.01 (Valid)	100.01
TCON-05-020	charges = 1999.99	TCOV-05-020	charges = 1999.99 (Valid)	1999.99
TCON-05-021	charges = 2000.00	TCOV-05-021	charges = 2000 (Valid)	2000
TCON-05-022	charges = 2000.01	TCOV-05-022	charges = 2000.01 (Invalid)	2000.01

Project Title: KHS Independent	Verification and Validation
--------------------------------	-----------------------------



TCON-05-004	Valid roomNo	TCOV-05-004	Valid roomNo (Valid output)	101
TCON-05-005	Invalid roomNo	TCOV-05-005	Invalid roomNo (Invalid output)	000
TCON-05-009	Valid charges	TCOV-05-009	Valid charges (Valid output)	900.0
TCON-05-010	Invalid charges	TCOV-05-010	Invalid charges (Invalid output)	2500.0

Table 2.2.5.2b Input Boundary Value Analysis Test Condition & Coverage

Test Plan ID: KHS\_TDS\_1\_1.0.0



## 2.2.5.3 Decision Table Testing

Test Condition ID	Test Condition
TCON-05-023	roomNo
TCON-05-024	charges
TCON-05-025	Display "Invalid roomNo"
TCON-05-026	Display "Invalid charges"
TCON-05-027	Display "Room successfully added"

Table 2.2.5.3a Add Room Decision Table Test Condition

Condition	TCOV-05-00 4	TCOV-05-00 9	TCOV-05-02 6	TCOV-05-02 7	TCOV-05-02 8
roomNo	roomNo < 1	1 ≤ roomNo ≤ 399	400 ≤ roomNo ≤ 500	1 ≤ roomNo ≤ 399	1 ≤ roomNo ≤ 399
charges	-	0 ≤ charges ≤ 99.99	-	2000.01 ≤ charges ≤ 3000	100≤ charges ≤ 2000
Output					
Display "Invalid roomNo"	×		×		
Display "Invalid charges"		*		×	
Display "Room successfully added"					×

Project Title: KHS Independent Verification and Validation			TERON O
	Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	LOS SOL

## **Table 2.2.5.3b Add Room Decision Table**

## 2.2.5.4 Use Case Testing

Use Case ID	UC005		
Use Case	F005 Add Room		
Purpose	To add	d room for database	
Actor	Hospi	tal Administrator	
Description	The hospital administrator starts by providing the information needed to add room into the database. If all valid requirements are input, the room will be added successfully into the database.		
Precondition	Syste	m prompt options	
Post-condition	Room	successfully added into database	
Scenario Name	Step	Action	
Main Flow	1	User choose option '1'	
	2	User enter roomNo	
	3	User enter charges	
	4	System validates roomNo	
	5	System validate charges	
	6	System add room into database	
Alternate Flow -	4.1	User Input invalid roomNo	
Invalid roomNo	4.2	System display error messages "Invalid room number"	
	4.3	The use case end	
Alternate Flow -	5.1	User Input invalid charges	
Invalid charges	5.2	System display error messages "Invalid charges"	
	5.3	The use case end	

Project Title: KHS Independent	STERON D	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	SINCE 202

Rules	None
Note	None

Table 2.2.5.4a Add Room Use Case



Test Condition ID	Test Condition	Test Coverage ID	Test Coverage	Test Data
TCON-05-028	Main Flow	TCOV-05-029	Main Flow	Valid roomNo = 123 Valid charges = 300
TCON-05-029	Alternate Flow - Invalid roomNo	TCOV-05-030	Alternate Flow - Invalid roomNo	Valid roomNo = 223 Valid charges = 300
TCON-05-030	Alternate Flow - Invalid charges	TCOV-05-031	Alternate Flow - Invalid charges	Valid roomNo = 123 Valid charges = 55

Table 2.2.5.4b Session Verification Use Case Test Conditions & Coverage

Project Title: KHS Independent	MERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	AND CE SON

Project Title: KHS Independent	MERON O	
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	AND CE SON

Project Title: KHS Independent Verification and Validation		
Date:06/04/2024	Test Plan ID:	

Test Plan ID: KHS\_TDS\_1\_1.0.0



## 2.3 Test Identification

### 2.3.1 Add Doctor

Test Case ID	Description	Test Coverage ID	Test ConditionID
TC-01-001	doctorId = 4999	TCOV-01-001	TCON-01-001
	gender = F	TCOV-01-005	TCON-01-005
		TCOV-01-010	TCON-01-010
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-023	TCON-01-023
TC-01-002	doctorId = 4999	TCOV-01-001	TCON-01-001
	gender = M	TCOV-01-005	TCON-01-005
		TCOV-01-010	TCON-01-010
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-023	TCON-01-023
TC-01-003	doctorId = 4999 gender = G	TCOV-01-001	TCON-01-001
		TCOV-01-005	TCON-01-005
		TCOV-01-007	TCON-01-007
		TCOV-01-006	TCON-01-006
		TCOV-01-009	TCON-01-009
		TCOV-01-023	TCON-01-023
		TCOV-01-024	TCON-01-024
TC-01-004	doctorId = 5000 gender = F	TCOV-01-002	TCON-01-002

Project Title: KHS Independent Verifi	ication and Validation
---------------------------------------	------------------------



		TCOV-01-004	TCON-01-004
		TCOV-01-011	TCON-01-011
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-007
		TCOV-01-022	TCON-01-022
		TCOV-01-016	TCON-01-016
		TCOV-01-017	TCON-01-017
		TCOV-01-021	TCON-01-021
TC-01-005	doctorId = 5000	TCOV-01-002	TCON-01-002
	gender = M	TCOV-01-004	TCON-01-004
		TCOV-01-011	TCON-01-011
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-022	TCON-01-022
		TCOV-01-016	TCON-01-016
		TCOV-01-018	TCON-01-018
		TCOV-01-021	TCON-01-021
TC-01-006	doctorId = 5000	TCOV-01-002	TCON-01-002
	gender = G	TCOV-01-004	TCON-01-004
		TCOV-01-011	TCON-01-011
		TCOV-01-006	TCON-01-006
		TCOV-01-024	TCON-01-024
		TCOV-01-016	TCON-01-016
		TCOV-01-020	TCON-01-020



		1	1
TC-01-007	doctorId = 5001 gender = F	TCOV-01-002	TCON-01-002
		TCOV-01-004	TCON-01-004
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-012	TCON-01-012
		TCOV-01-016	TCON-01-016
		TCOV-01-017	TCON-01-017
		TCOV-01-021	TCON-01-021
		TCOV-01-022	TCON-01-022
TC-01-008	doctorId = 5001	TCOV-01-002	TCON-01-002
	gender = M	TCOV-01-004	TCON-01-004
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-012	TCON-01-012
		TCOV-01-016	TCON-01-016
		TCOV-01-018	TCON-01-018
		TCOV-01-021	TCON-01-021
		TCOV-01-022	TCON-01-022
	doctorId = 5001 gender = G	TCOV-01-002	TCON-01-002
		TCOV-01-004	TCON-01-004
		TCOV-01-006	TCON-01-006
		TCOV-01-009	TCON-01-009
		TCOV-01-012	TCON-01-012
		TCOV-01-016	TCON-01-016



TCOV-01-020 TCON-01-020 TCOV-01-024 TCON-01-024  TC-01-010 doctorId = 7998 gender = F  TCOV-01-002 TCON-01-002 TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-017 TCOV-01-021 TCON-01-021 TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-008 TCON-01-002 TCOV-01-009 TCON-01-002 TCOV-01-009 TCON-01-002 TCOV-01-009 TCON-01-009 TCOV-01-009 TCON-01-009 TCOV-01-009 TCON-01-009 TCOV-01-013 TCON-01-013 TCOV-01-013 TCON-01-018 TCOV-01-018 TCON-01-018 TCOV-01-018 TCON-01-018				
TC-01-010  doctorId = 7998 gender = F  TCOV-01-002 TCOV-01-004 TCON-01-004 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCOV-01-016 TCOV-01-017 TCOV-01-017 TCOV-01-021 TCOV-01-021 TCOV-01-022 TCON-01-022 TC-01-011  doctorId = 7998 gender = M  TCOV-01-002 TCOV-01-003 TCOV-01-004 TCOV-01-005 TCOV-01-007 TCOV-01-007 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-013 TCOV-01-013 TCOV-01-016 TCOV-01-018 TCON-01-018			TCOV-01-020	TCON-01-020
TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-017 TCON-01-017 TCOV-01-021 TCON-01-021 TCOV-01-022 TCON-01-022 TC-01-011 doctorId = 7998 gender = M TCOV-01-004 TCON-01-004 TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-016 TCON-01-018			TCOV-01-024	TCON-01-024
TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-017 TCON-01-017  TCOV-01-021 TCON-01-021  TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-004 TCON-01-002  TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-016 TCON-01-018	TC-01-010		TCOV-01-002	TCON-01-002
TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-017 TCON-01-017  TCOV-01-021 TCON-01-021  TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-002 TCON-01-002  TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-018 TCON-01-018		gender = F	TCOV-01-004	TCON-01-004
TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-017 TCON-01-017  TCOV-01-021 TCON-01-021  TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-002 TCON-01-002  TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-018 TCON-01-018			TCOV-01-007	TCON-01-007
TCOV-01-016 TCON-01-016 TCOV-01-017 TCON-01-017 TCOV-01-021 TCON-01-021 TCOV-01-022 TCON-01-022 TC-01-011 doctorId = 7998 gender = M TCOV-01-002 TCON-01-002 TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-008	TCON-01-008
TCOV-01-017 TCON-01-017 TCOV-01-021 TCON-01-021 TCOV-01-022 TCON-01-022 TC-01-011 doctorId = 7998 gender = M TCOV-01-002 TCON-01-002 TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-013	TCON-01-013
TCOV-01-021 TCON-01-021 TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-002 TCON-01-002 TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-016	TCON-01-016
TCOV-01-022 TCON-01-022  TC-01-011 doctorId = 7998 gender = M  TCOV-01-002 TCON-01-002  TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-018 TCON-01-018			TCOV-01-017	TCON-01-017
TC-01-011			TCOV-01-021	TCON-01-021
gender = M  TCOV-01-004 TCON-01-004  TCOV-01-007 TCON-01-007  TCOV-01-008 TCON-01-008  TCOV-01-013 TCON-01-013  TCOV-01-016 TCON-01-016  TCOV-01-018 TCON-01-018			TCOV-01-022	TCON-01-022
TCOV-01-004 TCON-01-004 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018	TC-01-011		TCOV-01-002	TCON-01-002
TCOV-01-008 TCON-01-008 TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-004	TCON-01-004
TCOV-01-013 TCON-01-013 TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-007	TCON-01-007
TCOV-01-016 TCON-01-016 TCOV-01-018 TCON-01-018			TCOV-01-008	TCON-01-008
TCOV-01-018 TCON-01-018			TCOV-01-013	TCON-01-013
			TCOV-01-016	TCON-01-016
TCOV-01-021 TCON-01-021			TCOV-01-018	TCON-01-018
·			TCOV-01-021	TCON-01-021
TCOV-01-022 TCON-01-022			TCOV-01-022	TCON-01-022
TC-01-012 doctorId = 7998 TCOV-01-002 TCON-01-002	TC-01-012	doctorId = 7998 gender = G	TCOV-01-002	TCON-01-002
gender = G			TCOV-01-004	TCON-01-004
TCOV-01-006 TCON-01-006			TCOV-01-006	TCON-01-006
TCOV-01-009 TCON-01-009			TCOV-01-009	TCON-01-009



		TCOV-01-013	TCON-01-013
	TCOV-01-016	TCON-01-016	
		TCOV-01-020	TCON-01-020
		TCOV-01-024	TCON-01-024
TC-01-013	doctorId = 7999 gender = F	TCOV-01-003	TCON-01-003
		TCOV-01-004	TCON-01-004
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-014	TCON-01-014
		TCOV-01-016	TCON-01-016
		TCOV-01-017	TCON-01-017
		TCOV-01-021	TCON-01-021
		TCOV-01-022	TCON-01-022
TC-01-014	doctorld = 7999	TCOV-01-003	TCON-01-003
	gender = M	TCOV-01-004	TCON-01-004
		TCOV-01-007	TCON-01-007
		TCOV-01-008	TCON-01-008
		TCOV-01-014	TCON-01-014
		TCOV-01-016	TCON-01-016
		TCOV-01-018	TCON-01-018
		TCOV-01-021	TCON-01-021
		TCOV-01-022	TCON-01-022
TC-01-015	doctorld = 7999	TCOV-01-002	TCON-01-002
	gender = G	TCOV-01-004	TCON-01-004

Project Title: KHS Independent Verifi	ication and Validation
---------------------------------------	------------------------



TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-014 TCON-01-014 TCOV-01-016 TCON-01-016 TCOV-01-020 TCON-01-020 TCOV-01-024 TCON-01-024 TCOV-01-024 TCON-01-024 TCOV-01-024 TCON-01-024 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-022 TCOV-01-005 TCON-01-005 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-005 TCOV-01-008 TCON-01-005 TCOV-01-009 TCON-01-003 TCOV-01-005 TCON-01-005				
TCOV-01-014 TCON-01-014 TCOV-01-016 TCON-01-016 TCOV-01-020 TCON-01-020 TCOV-01-024 TCON-01-024 TCOV-01-024 TCON-01-024 TCOV-01-024 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-015 TCON-01-008 TCOV-01-023 TCON-01-022 TCOV-01-023 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-008 TCON-01-008 TCOV-01-008 TCON-01-008 TCOV-01-008 TCON-01-008 TCOV-01-009 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-009 TCON-01-009			TCOV-01-006	TCON-01-006
TCOV-01-016 TCON-01-016 TCOV-01-020 TCON-01-020 TCOV-01-024 TCON-01-024 TCOV-01-024 TCON-01-024 TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-003 TCOV-01-023 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-008 TCOV-01-023 TCON-01-008 TCOV-01-005 TCON-01-008 TCOV-01-006 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-009		TCOV-01-009	TCON-01-009	
TCOV-01-020 TCON-01-020 TCOV-01-024 TCON-01-024  TC-01-016 doctorId = 8001 gender = F  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-017 doctorId = 8001 gender = M  TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-008 TCOV-01-023 TCON-01-008 TCOV-01-015 TCON-01-005 TCOV-01-006 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-009		TCOV-01-014	TCON-01-014	
TC-01-016    Cov-01-024   TCON-01-024     TCOV-01-024   TCON-01-024     TCOV-01-003   TCON-01-003     TCOV-01-005   TCON-01-005     TCOV-01-007   TCON-01-007     TCOV-01-008   TCON-01-008     TCOV-01-015   TCON-01-015     TCOV-01-023   TCON-01-022     TCOV-01-007   TCON-01-003     TCOV-01-008   TCON-01-003     TCOV-01-005   TCON-01-005     TCOV-01-007   TCON-01-007     TCOV-01-008   TCON-01-007     TCOV-01-008   TCON-01-008     TCOV-01-015   TCON-01-015     TCOV-01-023   TCON-01-005     TCOV-01-023   TCON-01-005     TCOV-01-023   TCON-01-005     TCOV-01-006   TCON-01-005     TCOV-01-006   TCON-01-006     TCOV-01-009   TCON-01-009     TCOV-01-009   TCON-01-009     TCOV-01-0015   TCON-01-0015     TCOV-01			TCOV-01-016	TCON-01-016
TC-01-016  doctorId = 8001 gender = F  TCOV-01-003 TCON-01-003 TCON-01-005 TCON-01-005 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCOV-01-015 TCOV-01-023 TCON-01-022  TC-01-017  doctorId = 8001 gender = M  TCOV-01-003 TCOV-01-003 TCON-01-003 TCON-01-003 TCON-01-005 TCOV-01-005 TCOV-01-005 TCOV-01-006 TCOV-01-015 TCOV-01-023 TCON-01-008 TCON-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-008 TCOV-01-009 TCOV-01-006 TCOV-01-006 TCOV-01-009 TCOV-01-009 TCOV-01-009			TCOV-01-020	TCON-01-020
Gender = F			TCOV-01-024	TCON-01-024
TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-017 doctorId = 8001 gender = M  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-015 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-015 TCOV-01-023 TCON-01-002  TC-01-018 doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-009	TC-01-016		TCOV-01-003	TCON-01-003
TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-017 doctorId = 8001 gender = M  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-018 doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-005 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-015		gender = F	TCOV-01-005	TCON-01-005
TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-017 doctorId = 8001 gender = M  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022  TC-01-018 doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-015			TCOV-01-007	TCON-01-007
TCOV-01-023 TCON-01-022  TC-01-017  doctorId = 8001 gender = M  TCOV-01-003 TCON-01-003  TCOV-01-005 TCON-01-005  TCOV-01-007 TCON-01-007  TCOV-01-015 TCON-01-015  TCOV-01-023 TCON-01-015  TCOV-01-023 TCON-01-022  TC-01-018  doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003  TCOV-01-005 TCON-01-003  TCOV-01-005 TCON-01-005  TCOV-01-006 TCON-01-006  TCOV-01-009 TCON-01-009  TCOV-01-015 TCON-01-015			TCOV-01-008	TCON-01-008
TC-01-017  doctorId = 8001 gender = M  TCOV-01-003 TCON-01-003 TCON-01-005 TCOV-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCOV-01-015 TCOV-01-015 TCOV-01-023 TCON-01-022  TC-01-018  doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003 TCON-01-003 TCON-01-003 TCON-01-005 TCOV-01-006 TCOV-01-009 TCOV-01-009 TCOV-01-015 TCON-01-009			TCOV-01-015	TCON-01-015
TCOV-01-005   TCON-01-005     TCOV-01-007   TCON-01-007     TCOV-01-008   TCON-01-008     TCOV-01-015   TCON-01-015     TCOV-01-023   TCON-01-022     TC-01-018   doctorId = 8001 gender = G     TCOV-01-003   TCON-01-003     TCOV-01-005   TCON-01-005     TCOV-01-006   TCON-01-006     TCOV-01-009   TCON-01-009     TCOV-01-015   TCON-01-015			TCOV-01-023	TCON-01-022
TCOV-01-005 TCON-01-005 TCOV-01-007 TCON-01-007 TCOV-01-008 TCON-01-008 TCOV-01-015 TCON-01-015 TCOV-01-023 TCON-01-022 TC-01-018 doctorId = 8001 gender = G TCOV-01-003 TCON-01-003 TCOV-01-005 TCON-01-005 TCOV-01-006 TCON-01-006 TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-015	TC-01-017		TCOV-01-003	TCON-01-003
TCOV-01-008 TCON-01-008  TCOV-01-015 TCON-01-015  TCOV-01-023 TCON-01-022  TC-01-018 doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003  TCOV-01-005 TCON-01-005  TCOV-01-006 TCON-01-006  TCOV-01-009 TCON-01-009  TCOV-01-015 TCON-01-015		gender = M	TCOV-01-005	TCON-01-005
TCOV-01-015 TCON-01-015  TCOV-01-023 TCON-01-022  TC-01-018 doctorId = 8001 gender = G  TCOV-01-003 TCON-01-003 TCON-01-005  TCOV-01-005 TCON-01-005  TCOV-01-006 TCON-01-009  TCOV-01-015 TCON-01-015			TCOV-01-007	TCON-01-007
			TCOV-01-008	TCON-01-008
$ \begin{array}{c ccccccccccccccccccccccccccccccccccc$			TCOV-01-015	TCON-01-015
gender = G  TCOV-01-005 TCON-01-005  TCOV-01-006 TCON-01-006  TCOV-01-009 TCON-01-009  TCOV-01-015 TCON-01-015			TCOV-01-023	TCON-01-022
TCOV-01-005 TCON-01-005  TCOV-01-006 TCON-01-006  TCOV-01-009 TCON-01-009  TCOV-01-015 TCON-01-015	TC-01-018		TCOV-01-003	TCON-01-003
TCOV-01-009 TCON-01-009 TCOV-01-015 TCON-01-015		gender = G	TCOV-01-005	TCON-01-005
TCOV-01-015 TCON-01-015			TCOV-01-006	TCON-01-006
			TCOV-01-009	TCON-01-009
TCOV-01-023 TCON-01-023			TCOV-01-015	TCON-01-015
			TCOV-01-023	TCON-01-023



		TCOV-01-024	TCON-01-024
TC-01-019	doctorId = 9	TCOV-01-016	TCON-01-016
	gender = Don't care	TCOV-01-017	TCON-01-017
		TCOV-01-018	TCON-01-018
		TCOV-01-019	TCON-01-019
		TCOV-01-023	TCON-01-023
TC-01-020	doctorld = 10005	TCOV-01-016	TCON-01-016
	gender = Don't care	TCOV-01-017	TCON-01-017
		TCOV-01-018	TCON-01-018
		TCOV-01-019	TCON-01-019
		TCOV-01-023	TCON-01-023
TC-01-021	Add Doctor Main Flow	TCOV-01-002	TCON-01-018
		TCOV-01-004	TCON-01-018
		TCOV-01-007	TCON-01-018
		TCOV-01-008	TCON-01-018
		TCOV-01-012	TCON-01-018
		TCOV-01-013	TCON-01-018
		TCOV-01-014	TCON-01-018
		TCOV-01-016	TCON-01-018
		TCOV-01-017	TCON-01-018
		TCOV-01-018	TCON-01-018
		TCOV-01-021	TCON-01-018
		TCOV-01-022	TCON-01-018
TC-01-022	Add Doctor	TCOV-01-010	TCON-01-010

Project Title: KHS Independent Verification and Validation		STERON O
Date:06/04/2024	Test Plan ID: KHS_TDS_1_1.0.0	INCE SOL

	AlternateFlow -	TCOV-01-011	TCON-01-011
	Invalid doctorld	TCOV-01-007	TCON-01-007
		TCOV-01-009	TCON-01-009
		TCOV-01-008	TCON-01-008
		TCOV-01-002	TCON-01-002
		TCOV-01-015	TCON-01-015
TC-01-023	Add Doctor	TCOV-01-011	TCON-01-011
	AlternateFlow - Invalid gender	TCOV-01-004	TCON-01-004
		TCOV-01-006	TCON-01-006
		TCOV-01-009	TCON-01-009
		TCOV-01-016	TCON-01-016
		TCOV-01-019	TCON-01-019
		TCOV-01-024	TCON-01-024

**Table 2.3.1 Add Doctor Test Identification** 

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.3.2 Add Patient

Test Case ID	Description	Test Coverage ID	Test Condition ID
TC-02-001	Patient ID = 999	TCOV-02-002	TCON-02-002
	IllnessLength = 0	TCOV-02-005	TCON-02-005
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-009
		TCOV-02-010	TCON-02-010
		TCOV-02-018	TCON-02-018
TC-02-002	Patient ID = 999	TCOV-02-002	TCON-02-002
	IllnessLength = 1	TCOV-02-005	TCON-02-005
		TCOV-02-006	TCON-02-006
		TCOV-02-008	TCON-02-008
		TCOV-02-010	TCON-02-010
		TCOV-02-019	TCON-02-019
TC-02-003		TCOV-02-001	TCON-02-001
	IllnessLength = 0	TCOV-02-004	TCON-02-004
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-00
		TCOV-02-011	TCON-02-011
		TCOV-02-018	TCON-02-018
TC-02-004	Patient ID = 1000	TCOV-02-001	TCON-02-001
	IllnessLength = 1	TCOV-02-004	TCON-02-004
		TCOV-02-006	TCON-02-006

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-02-008	TCON-02-008
		TCOV-02-011	TCON-02-011
		TCOV-02-019	TCON-02-019
TC-02-005	Patient ID = 1001	TCOV-02-002	TCON-02-002
	IllnessLength = 0	TCOV-02-004	TCON-02-004
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-009
		TCOV-02-012	TCON-02-012
		TCOV-02-018	TCON-02-018
TC-02-006	Patient ID = 1001	TCOV-02-001	TCON-02-001
	IllnessLength = 1	TCOV-02-004	TCON-02-004
		TCOV-02-006	TCON-02-006
		TCOV-02-008	TCON-02-008
		TCOV-02-012	TCON-02-012
		TCOV-02-019	TCON-02-019
TC-02-007	Patient ID = 3998	TCOV-02-002	TCON-02-002
	IllnessLength = 0	TCOV-02-005	TCON-02-005
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-009
		TCOV-02-013	TCON-02-013
		TCOV-02-018	TCON-02-018
TC-02-008	Patient ID = 3998	TCOV-02-002	TCON-02-002
IllnessLe	IllnessLength = 1	TCOV-02-005	TCON-02-005
		TCOV-02-006	TCON-02-006

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-02-008	TCON-02-008
		TCOV-02-013	TCON-02-013
		TCOV-02-019	TCON-02-019
TC-02-009	Patient ID = 3999	TCOV-02-001	TCON-02-001
	IllnessLength = 0	TCOV-02-004	TCON-02-004
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-00
		TCOV-02-014	TCON-02-014
		TCOV-02-018	TCON-02-018
TC-02-010	Patient ID = 3999	TCOV-02-001	TCON-02-001
	IllnessLength = 1	TCOV-02-004	TCON-02-004
		TCOV-02-006	TCON-02-006
		TCOV-02-008	TCON-02-008
		TCOV-02-014	TCON-02-014
		TCOV-02-019	TCON-02-019
TC-02-011	Patient ID = 4000	TCOV-02-002	TCON-02-002
	IllnessLength = 0	TCOV-02-004	TCON-02-004
		TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-009
		TCOV-02-015	TCON-02-015
		TCOV-02-018	TCON-02-018
TC-02-012	Patient ID = 4000	TCOV-02-001	TCON-02-001
	IllnessLength = 1	TCOV-02-004	TCON-02-004
		TCOV-02-006	TCON-02-006
	<del></del>		

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



		TCOV-02-008	TCON-02-008
		TCOV-02-015	TCON-02-015
		TCOV-02-019	TCON-02-019
TC-02-013	Patient ID = 560	TCOV-02-022	TCON-02-022
	IllnessLength = don't care	TCOV-02-023	TCON-02-023
		TCOV-02-024	TCON-02-024
TC-02-014	Patient ID = 5678	TCOV-02-022	TCON-02-022
	IllnessLength = don't care	TCOV-02-023	TCON-02-023
		TCOV-02-024	TCON-02-024
TC-02-015	Main Flow	TCOV-02-001	TCON-02-001
		TCOV-02-004	TCON-02-004
		TCOV-02-006	TCON-02-006
		TCOV-02-008	TCON-02-008
		TCOV-02-027	TCON-02-027
TC-02-016	Alternate Flow -	TCOV-02-002	TCON-02-002
	Invalid patientID	TCOV-02-003	TCON-02-003
		TCOv-02-005	TCON-02-005
		TCOV-02-028	TCON-02-028
TC-02-017	Alternate Flow - Invalid Illness	TCOV-02-007	TCON-02-007
		TCOV-02-009	TCON-02-009
		TCOV-02-029	TCON-02-029

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.3.3 Add Medicine

Test Case ID	Description	Test Coverage ID	Test Condition ID
TC-03-001	medicineCode = 3999	TCOV-03-002	TCON-03-002
	price = 0.00	TCOV-03-005	TCON-03-005
		TCOV-03-008	TCON-03-008
		TCOV-03-010	TCON-03-010
		TCOV-03-011	TCON-03-011
		TCOV-03-017	TCON-03-017
TC-03-002	medicineCode = 3999	TCOV-03-002	TCON-03-002
	price = 1.00	TCOV-03-004	TCON-03-004
		TCOV-03-008	TCON-03-008
		TCOV-03-009	TCON-03-009
		TCOV-03-011	TCON-03-011
		TCOV-03-018	TCON-03-018
TC-03-003		TCOV-03-002	TCON-03-002
	price = 2.00	TCOV-03-004	TCON-03-004
		TCOV-03-008	TCON-03-008
		TCOV-03-009	TCON-03-009
		TCOV-03-011	TCON-03-011
		TCOV-03-019	TCON-03-019
TC-03-004	medicineCode = 3999	TCOV-03-002	TCON-03-002
	price = 49.00	TCOV-03-004	TCON-03-004
		TCOV-03-008	TCON-03-008

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-03-009	TCON-03-009
		TCOV-03-011	TCON-03-019
		TCOV-03-020	TCON-03-020
TC-03-005	medicineCode = 3999	TCOV-03-002	TCON-03-002
	price = 50.00	TCOV-03-004	TCON-03-004
		TCOV-03-008	TCON-03-008
		TCOV-03-009	TCON-03-009
		TCOV-03-011	TCON-03-011
		TCOV-03-021	TCON-03-021
TC-03-006	medicineCode = 3999	TCOV-03-002	TCON-03-002
	price = 51.00	TCOV-03-005	TCON-03-005
		TCOV-03-008	TCON-03-008
		TCOV-03-010	TCON-03-010
		TCOV-03-011	TCON-03-011
		TCOV-03-022	TCON-03-022
TC-03-007	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price = 0.00	TCOV-03-005	TCON-03-005
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-012	TCON-03-012
		TCOV-03-017	TCON-03-017
TC-03-008	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price = 1.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-03-009	TCON-03-009
		TCOV-03-012	TCON-03-012
		TCOV-03-018	TCON-03-018
TC-03-009	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price = 2.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-012	TCON-03-012
		TCOV-03-019	TCON-03-019
TC-03-010	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price =49.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-012	TCON-03-012
		TCOV-03-020	TCON-03-020
TC-03-011	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price = 50.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-012	TCON-03-012
		TCOV-03-021	TCON-03-021
TC-03-012	medicineCode = 4000	TCOV-03-001	TCON-03-001
	price = 51.00	TCOV-03-006	TCON-03-006
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-03-010	TCON-03-010
		TCOV-03-012	TCON-03-012
		TCOV-03-022	TCON-03-022
TC-03-013	medicine Code = 4001	TCOV-03-001	TCON-03-001
	price = 0.00	TCOV-03-005	TCON-03-005
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-013	TCON-03-013
		TCOV-03-017	TCON-03-017
TC-03-014	medicine Code = 4001	TCOV-03-001	TCON-03-001
	price = 1.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-013	TCON-03-013
		TCOV-03-018	TCON-03-018
TC-03-015	medicine Code = 4001	TCOV-03-001	TCON-03-001
	price = 2.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-013	TCON-03-013
		TCOV-03-019	TCON-03-019
TC-03-016	medicine Code = 4001	TCOV-03-001	TCON-03-001
	price =49.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

WINCE 202

		TCOV-03-009	TCON-03-009
		TCOV-03-013	TCON-03-013
		TCOV-03-020	TCON-03-020
TC-03-017	medicine Code = 4001	TCOV-03-001	TCON-03-001
	price = 50.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-013	TCON-03-013
		TCOV-03-021	TCON-03-021
TC-03-018	medicine Code = 4001 price = 51.00	TCOV-03-001	TCON-03-001
		TCOV-03-006	TCON-03-006
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-013	TCON-03-013
		TCOV-03-022	TCON-03-022
TC-03-019	medicineCode = 4513	TCOV-03-001	TCON-03-001
	price = 0.00	TCOV-03-005	TCON-03-005
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-014	TCON-03-014
		TCOV-03-017	TCON-03-017
TC-03-020	medicineCode = 4513 price = 1.00	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



		TCOV-03-009	TCON-03-009
		TCOV-03-014	TCON-03-014
		TCOV-03-018	TCON-03-018
TC-03-021	medicineCode = 4513 price = 2.00	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-014	TCON-03-014
		TCOV-03-019	TCON-03-019
TC-03-022	medicineCode = 4513 price =49.00	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-014	TCON-03-014
		TCOV-03-020	TCON-03-020
TC-03-023	medicineCode = 4513	TCOV-03-001	TCON-03-001
	price = 50.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-014	TCON-03-014
		TCOV-03-021	TCON-03-021
TC-03-024	medicineCode = 4513 price = 51.00	TCOV-03-001	TCON-03-001
		TCOV-03-006	TCON-03-006
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0



		TCOV-03-010	TCON-03-010
		TCOV-03-014	TCON-03-014
		TCOV-03-022	TCON-03-022
TC-03-025	medicineCode = 4514	TCOV-03-001	TCON-03-001
	price = 0.00	TCOV-03-005	TCON-03-005
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-015	TCON-03-015
		TCOV-03-017	TCON-03-017
TC-03-026	medicineCode = 4514 price = 1.00	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-015	TCON-03-015
		TCOV-03-018	TCON-03-018
TC-03-027	medicineCode = 4514	TCOV-03-001	TCON-03-001
	price = 2.00	TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-015	TCON-03-015
		TCOV-03-019	TCON-03-019
TC-03-028	medicineCode = 4514 price =49.00	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



TCOV-03-009 TCON-03-00  TCOV-03-015 TCON-03-01  TCOV-03-020 TCON-03-02  TCOV-03-020 TCON-03-02  TCOV-03-001 TCON-03-00  TCOV-03-004 TCON-03-00  TCOV-03-007 TCON-03-00  TCOV-03-009 TCON-03-00  TCOV-03-015 TCON-03-01  TCOV-03-015 TCON-03-02	5 ) 1 1 7
TC-03-029 medicineCode = 4514 price = 50.00 TCOV-03-020 TCON-03-020 TCON-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-000 TCOV-03-015 TCON-03-01	) 1 4 7
TC-03-029 medicineCode = 4514 price = 50.00 TCOV-03-001 TCON-03-000 TCOV-03-007 TCON-03-000 TCOV-03-009 TCON-03-000 TCOV-03-015 TCON-03-010 TCOV-03-015 TCON-03-010 TCOV-03-015	1 7
TCOV-03-004 TCON-03-004 TCOV-03-007 TCON-03-000 TCOV-03-009 TCON-03-000 TCOV-03-015 TCON-03-019	1 7 9
TCOV-03-004 TCON-03-004 TCOV-03-007 TCON-03-000 TCOV-03-009 TCON-03-000 TCOV-03-015 TCON-03-01	7
TCOV-03-009 TCON-03-009 TCOV-03-015 TCON-03-019	9
TCOV-03-015 TCON-03-01	
	5
TCOV-03-021 TCON-03-02	
	l
TC-03-030 medicineCode = 4514 TCOV-03-001 TCON-03-00	l
price = 51.00 TCOV-03-006 TCON-03-00	3
TCOV-03-007 TCON-03-00	7
TCOV-03-010 TCON-03-010	)
TCOV-03-015 TCON-03-01	5
TCOV-03-022 TCON-03-02	2
TC-03-031 medicineCode = 4515 TCOV-03-003 TCON-03-00	3
price = 0.00 TCOV-03-005 TCON-03-00	5
TCOV-03-008 TCON-03-00	3
TCOV-03-010 TCON-03-010	)
TCOV-03-016 TCON-03-01	3
TCOV-03-017 TCON-03-01	7
TC-03-032 medicineCode = 4515 TCOV-03-003 TCON-03-00	3
price = 1.00 TCOV-03-004 TCON-03-004	1
TCOV-03-008 TCON-03-00	

Date:06/04/2024 Test Plan ID:



TCOV-03-009 TCON-03-009  TCOV-03-016 TCON-03-016  TCOV-03-018 TCON-03-018  TCOV-03-018 TCON-03-018  TCOV-03-003 TCON-03-003  TCOV-03-004 TCON-03-004  TCOV-03-008 TCON-03-008  TCOV-03-009 TCON-03-009  TCOV-03-016 TCON-03-016  TCOV-03-019 TCON-03-019  TCOV-03-034 medicineCode = 4515 TCOV-03-003  TCOV-03-003 TCON-03-003
TC-03-033 medicineCode = 4515 price = 2.00 TCOV-03-003 TCON-03-003 TCON-03-004 TCON-03-004 TCON-03-008 TCOV-03-008 TCOV-03-008 TCOV-03-008 TCOV-03-008 TCOV-03-008 TCOV-03-016 TCOV-03-016 TCOV-03-018
TC-03-033 medicineCode = 4515 price = 2.00 TCOV-03-003 TCON-03-003 TCON-03-004 TCON-03-004 TCON-03-008 TCOV-03-008 TCOV-03-009 TCON-03-008 TCOV-03-016 TCOV-03-016 TCOV-03-019 TCON-03-019
TCOV-03-004 TCON-03-004 TCOV-03-008 TCON-03-008 TCOV-03-009 TCON-03-008 TCOV-03-016 TCON-03-018
TCOV-03-004 TCON-03-004 TCOV-03-008 TCON-03-008 TCOV-03-009 TCON-03-009 TCOV-03-016 TCON-03-019 TCOV-03-019 TCON-03-019
TCOV-03-009 TCON-03-009 TCOV-03-016 TCON-03-019 TCOV-03-019 TCON-03-019
TCOV-03-016 TCON-03-016 TCOV-03-019 TCON-03-019
TCOV-03-019 TCON-03-019
TC-03-034 medicineCode = 4515 TCOV-03-003 TCON-03-003
price =49.00 TCOV-03-004 TCON-03-004
TCOV-03-008 TCON-03-008
TCOV-03-009 TCON-03-009
TCOV-03-016 TCON-03-016
TCOV-03-020 TCON-03-020
TC-03-035 medicineCode = 4515 TCOV-03-003 TCON-03-003
price = 50.00 TCOV-03-004 TCON-03-004
TCOV-03-008 TCON-03-008
TCOV-03-009 TCON-03-009
TCOV-03-016 TCON-03-016
TCOV-03-021 TCON-03-02 <sup>-2</sup>
TC-03-036 medicineCode = 4515 TCOV-03-003 TCON-03-003
price = 51.00 TCOV-03-006 TCON-03-006
TCOV-03-008 TCON-03-008

Date:06/04/2024 Test Plan ID:



		TCOV-03-010	TCON-03-010
		TCOV-03-016	TCON-03-016
		TCOV-03-022	TCON-03-022
TC-03-037	medicineCode < 4000	TCOV-03-023	TCON-03-023
	price = "don't care"		TCON-03-025
TC-03-038	medicineCode > 4514	TCOV-03-027	TCON-03-023
	price = "don't care"		TCON-03-025
TC-03-039	Add Medicine Main Flow	TCOV-03-001	TCON-03-001
		TCOV-03-004	TCON-03-004
		TCOV-03-007	TCON-03-007
		TCOV-03-009	TCON-03-009
		TCOV-03-029	TCON-03-028
TC-03-040		TCOV-03-002	TCON-03-002
	Invalid Id	TCOV-03-003	TCON-03-003
		TCOV-03-008	TCON-03-008
		TCOV-03-030	TCON-03-029
TC-03-041	Add medicine AlternateFlow -	TCOV-03-001	TCON-03-001
	Invalid price	TCOV-03-005	TCON-03-005
		TCOV-03-006	TCON-03-006
		TCOV-03-007	TCON-03-007
		TCOV-03-010	TCON-03-010
		TCOV-03-031	TCON-03-030

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



# 2.3.4 Add Nurse

Test Case ID	Description	Test Coverage ID	Test Condition ID
TC-04-001	nurseld = 7999	TCOV-04-010	TCON-04-010
	gender = F	TCOV-04-004	TCON-04-004
		TCOV-04-007	TCON-04-007
		TCOV-04-008	TCON-04-008
		TCOV-04-023	TCON-04-02
		TCOV-04-002	TCON-04-002
TC-04-002	nurseld = 7999	TCOV-04-010	TCON-04-010
	gender = M	TCOV-04-004	TCON-04-004
		TCOV-04-007	TCON-04-007
		TCOV-04-008	TCON-04-008
		TCOV-04-023	TCON-04-023
		TCOV-04-002	TCON-04-002
TC-04-003		TCOV-04-010	TCON-04-010
	gender = G	TCOV-04-005	TCON-04-005
		TCOV-04-007	TCON-04-007
		TCOV-04-009	TCON-04-009
		TCOV-04-023	TCON-04-023
		TCOV-04-024	TCON-04-024
		TCOV-04-002	TCON-04-002
	nurseld = 8000	TCOV-04-011	TCON-04-011
	gender = F	TCOV-04-004	TCON-04-004

Date:06/04/2024 Test Plan ID

Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-018	TCON-04-018
		TCOV-04-021	TCON-04-021
TC-04-005	nurseld = 8000	TCOV-04-011	TCON-04-011
	gender = M	TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-017	TCON-04-017
		TCOV-04-021	TCON-04-021
TC-04-006	nurseld = 8000	TCOV-04-011	TCON-04-011
	gender = O	TCOV-04-005	TCON-04-005
		TCOV-04-006	TCON-04-006
		TCOV-04-009	TCON-04-009
		TCOV-04-024	TCON-04-024
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-020	TCON-04-020

Date:06/04/2024 Test Plan ID:



TC-04-007	nurseld = 8001	TCOV-04-012	TCON-04-012
	gender = F	TCOV-04-004	TCON-04-004
		TCOV-0-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-018	TCON-04-018
		TCOV-04-021	TCON-04-021
TC-04-008	nurseld = 8001	TCOV-04-012	TCON-04-012
	gender = M	TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-017	TCON-04-017
		TCOV-04-021	TCON-04-021
TC-04-009	nurseld = 8001	TCOV-04-012	TCON-04-012
	gender = A	TCOV-04-005	TCON-04-005
		TCOV-04-006	TCON-04-006
		TCOV-04-009	TCON-04-009
		TCOV-04-024	TCON-04-024
		TCOV-04-001	TCON-04-001

Date:06/04/2024 Test Plan ID:



		TCOV-04-016	TCON-04-016
		TCOV-04-020	TCON-04-020
TC-04-010	nurseld = 9998	TCOV-04-013	TCON-04-013
	gender = F	TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-018	TCON-04-018
		TCOV-04-021	TCON-04-021
TC-04-011	nurseld = 9998 gender = M	TCOV-04-013	TCON-04-013
		TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-017	TCON-04-017
		TCOV-04-021	TCON-04-021
TC-04-012	nurseld = 9998	TCOV-04-013	TCON-04-013
	gender = B	TCOV-04-005	TCON-04-005
		TCOV-04-006	TCON-04-006
		TCOV-04-009	TCON-04-009
		100 104 000	1001107000

Date:06/04/2024 Test Plan ID:



		TCOV-04-024	TCON-04-024
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-020	TCON-04-020
TC-04-013	nurseld = 9999	TCOV-04-014	TCON-04-014
	gender = F	TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-018	TCON-04-018
		TCOV-04-021	TCON-04-021
TC-04-014	nurseld = 9999	TCOV-04-014	TCON-04-014
	gender = M	TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
	TCOV-04-017	TCON-04-017	
		TCOV-04-021	TCON-04-021
TC-04-015	nurseld = 9999	TCOV-04-014	TCON-04-014
	gender = K	TCOV-04-005	TCON-04-005

Date:06/04/2024 Test Plan ID:



		TCOV-04-006	TCON-04-006
		TCOV-04-009	TCON-04-009
		TCOV-04-024	TCON-04-024
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-020	TCON-04-020
TC-04-016	nurseld = 10000	TCOV-04-015	TCON-04-015
	gender = F	TCOV-04-004	TCON-04-004
		TCOV-04-007	TCON-04-007
		TCOV-04-008	TCON-04-008
		TCOV-04-023	TCON-04-023
		TCOV-04-003	TCON-04-003
TC-04-017	nurseld = 10000 gender = M	TCOV-04-015	TCON-04-015
		TCOV-04-004	TCON-04-004
		TCOV-04-007	TCON-04-007
		TCOV-04-008	TCON-04-008
		TCOV-04-023	TCON-04-023
		TCOV-04-003	TCON-04-003
TC-04-018	nurseld = 10000	TCOV-04-015	TCON-04-015
	gender = Y	TCOV-04-005	TCON-04-005
		TCOV-04-007	TCON-04-007
		TCOV-04-009	TCON-04-009
		TCOV-04-023	TCON-04-023
		TCOV-04-024	TCON-04-024

Date:06/04/2024 Test Plan ID:



		TCOV-04-003	TCON-04-003
TC-04-019	nurseld = 9	TCOV-04-016	TCON-04-016
	gender = Don't care	TCOV-04-017	TCON-04-017
		TCOV-04-018	TCON-04-018
		TCOV-04-019	TCON-04-019
		TCOV-04-023	TCON-04-023
TC-04-020	nurseld = 10025	TCOV-04-016	TCON-04-016
	gender = Don't care	TCOV-04-017	TCON-04-017
		TCOV-04-018	TCON-04-019
		TCOV-04-019	TCON-04-019
		TCOV-04-023	TCON-04-023
TC-03-021	Add Nurse Main Flow	TCOV-04-011	TCON-04-011
		TCOV-04-004	TCON-04-004
		TCOV-04-006	TCON-04-006
		TCOV-04-008	TCON-04-008
		TCOV-04-022	TCON-04-022
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-018	TCON-04-018
		TCOV-04-021	TCON-04-021
		TCOV-04-017	TCON-04-017
		TCOV-04-012	TCON-04-012
		TCOV-04-013	TCON-04-013
		TCOV-04-014	TCON-04-014

Date:06/04/2024 Test Plan ID:



			_
TC-03-022	Add Nurse AlternateFlow -	TCOV-04-010	TCON-04-010
	Invalid nurseID	TCOV-04-004	TCON-04-004
		TCOV-04-007	TCON-04-007
		TCOV-04-008	TCON-04-008
		TCOV-04-023	TCON-04-02
		TCOV-04-002	TCON-04-002
		TCOV-04-015	TCON-04-015
TC-03-023	Add Nurse AlternateFlow -	TCOV-04-011	TCON-04-011
	Invalid gender	TCOV-04-005	TCON-04-005
		TCOV-04-006	TCON-04-006
		TCOV-04-009	TCON-04-009
		TCOV-04-024	TCON-04-024
		TCOV-04-001	TCON-04-001
		TCOV-04-016	TCON-04-016
		TCOV-04-020	TCON-04-020
		TCOV-04-012	TCON-04-012
		TCOV-04-013	TCON-04-013
		TCOV-04-014	TCON-04-014

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



# 2.3.5 Add Room

Test Case ID	Description	Test Coverage ID	Test Condition ID
TC-05-001	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 99.99	TCOV-05-005	TCON-05-005
		TCOV-05-006	TCON-05-006
		TCOV-05-010	TCON-05-010
		TCOV-05-011	TCON-05-011
		TCOV-05-017	TCON-05-017
TC-05-002	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 100.00	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-011	TCON-05-011
		TCOV-05-018	TCON-05-018
	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 100.01	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-011	TCON-05-011
		TCOV-05-019	TCON-05-019
TC-05-004	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 1999.99	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

SINCE SOL

		TCOV-05-009	TCON-05-009
		TCOV-05-011	TCON-05-011
		TCOV-05-020	TCON-05-020
TC-05-005	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 2000.00	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-011	TCON-05-011
		TCOV-05-021	TCON-05-021
TC-05-006	roomNo = 000	TCOV-05-001	TCON-05-001
	charges = 2000.01	TCOV-05-005	TCON-05-005
		TCOV-05-008	TCON-05-008
		TCOV-05-010	TCON-05-010
		TCOV-05-011	TCON-05-011
		TCOV-05-022	TCON-05-022
TC-05-007	roomNo = 001	TCOV-05-002	TCON-05-002
	charges = 99.99	TCOV-05-004	TCON-05-004
		TCOV-05-006	TCON-05-006
		TCOV-05-010	TCON-05-010
		TCOV-05-012	TCON-05-012
		TCOV-05-017	TCON-05-017
TC-05-008	TC-05-008 roomNo = 001 charges = 100.00	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-05-009	TCON-05-009
		TCOV-05-012	TCON-05-012
		TCOV-05-018	TCON-05-018
TC-05-009	roomNo = 001	TCOV-05-002	TCON-05-002
	charges = 100.01	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-012	TCON-05-012
		TCOV-05-019	TCON-05-019
TC-05-010	roomNo = 001	TCOV-05-002	TCON-05-002
	charges = 1999.99	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-012	TCON-05-012
		TCOV-05-020	TCON-05-020
TC-05-011	roomNo = 001 charges = 2000.00	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-012	TCON-05-012
		TCOV-05-021	TCON-05-021
TC-05-012	roomNo = 001	TCOV-05-002	TCON-05-002
	charges = 2000.01	TCOV-05-004	TCON-05-004
		TCOV-05-008	TCON-05-008

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-05-010	TCON-05-010
		TCOV-05-012	TCON-05-012
		TCOV-05-022	TCON-05-022
TC-05-013	roomNo = 002	TCOV-05-002	TCON-05-002
	charges = 99.99	TCOV-05-004	TCON-05-004
		TCOV-05-006	TCON-05-006
		TCOV-05-010	TCON-05-010
		TCOV-05-013	TCON-05-013
		TCOV-05-017	TCON-05-017
TC-05-014	roomNo = 002 charges = 100.00	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-013	TCON-05-013
		TCOV-05-018	TCON-05-018
TC-05-015	roomNo = 002 charges = 100.01	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-013	TCON-05-013
		TCOV-05-019	TCON-05-019
TC-05-016	roomNo = 002	TCOV-05-002	TCON-05-002
	charges = 1999.99	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-05-009	TCON-05-009
		TCOV-05-013	TCON-05-013
		TCOV-05-020	TCON-05-020
TC-05-017	roomNo = 002	TCOV-05-002	TCON-05-002
	charges = 2000.00	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-013	TCON-05-013
		TCOV-05-021	TCON-05-021
TC-05-018	roomNo = 002	TCOV-05-002	TCON-05-002
	charges = 2000.01	TCOV-05-004	TCON-05-004
		TCOV-05-008	TCON-05-008
		TCOV-05-010	TCON-05-010
		TCOV-05-013	TCON-05-013
		TCOV-05-022	TCON-05-022
TC-05-019	TC-05-019 roomNo = 398 charges = 99.99	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-006	TCON-05-006
		TCOV-05-010	TCON-05-010
		TCOV-05-014	TCON-05-014
		TCOV-05-017	TCON-05-017
TC-05-020	roomNo = 398	TCOV-05-002	TCON-05-002
	charges = 100.00	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID:



TCOV-05-009 TCON-05-009 TCOV-05-014 TCON-05-014 TCOV-05-018 TCON-05-018 TC-05-021 roomNo = 398 charges = 100.01 TCOV-05-004 TCON-05-002 TCOV-05-004 TCON-05-004 TCOV-05-009 TCON-05-009 TCOV-05-019 TCON-05-019 TCOV-05-019 TCON-05-019 TCOV-05-019 TCON-05-002 TCOV-05-004 TCON-05-002 TCOV-05-004 TCON-05-009 TCOV-05-019 TCON-05-009 TCOV-05-019 TCON-05-002 TCOV-05-004 TCON-05-002 TCOV-05-004 TCON-05-004 TCOV-05-007 TCON-05-007 TCOV-05-014 TCON-05-009 TCOV-05-014 TCON-05-014 TCOV-05-014 TCON-05-014 TCOV-05-020 TCON-05-020 TC-05-023 roomNo = 398 charges = 2000.00				
TC-05-021			TCOV-05-009	TCON-05-009
TC-05-021			TCOV-05-014	TCON-05-014
Charges = 100.01			TCOV-05-018	TCON-05-018
TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-009 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-019 TCON-05-019  TCOV-05-019 TCON-05-019  TCOV-05-019 TCON-05-019  TCOV-05-002 TCON-05-002  TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-009 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-020 TCON-05-020  TCOV-05-020 TCON-05-020  TCOV-05-020 TCON-05-020	TC-05-021		TCOV-05-002	TCON-05-002
TCOV-05-009 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-019 TCON-05-019  TCOV-05-019 TCON-05-019  TCOV-05-019 TCON-05-019  TCOV-05-002 TCON-05-002  TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-009 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-020 TCON-05-020  TC-05-023 roomNo = 398 Charges = 2000.00		charges = 100.01	TCOV-05-004	TCON-05-004
TCOV-05-014 TCON-05-014  TCOV-05-019 TCON-05-019  TC-05-022 roomNo = 398 charges = 1999.99  TCOV-05-004 TCON-05-002  TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-014 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-020 TCON-05-020  TCOV-05-020 TCON-05-020			TCOV-05-007	TCON-05-007
TC-05-022			TCOV-05-009	TCON-05-009
TC-05-022			TCOV-05-014	TCON-05-014
Charges = 1999.99  TCOV-05-004  TCOV-05-007  TCON-05-007  TCON-05-009  TCOV-05-014  TCOV-05-014  TCOV-05-020  TCON-05-020  TCON-05-020  TCON-05-020			TCOV-05-019	TCON-05-019
TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-009 TCON-05-009  TCOV-05-014 TCON-05-014  TCOV-05-020 TCON-05-020  TC-05-023 roomNo = 398 TCOV-05-002 TCON-05-002  charges = 2000.00	TC-05-022		TCOV-05-002	TCON-05-002
TCOV-05-009 TCON-05-009 TCOV-05-014 TCON-05-014 TCOV-05-020 TCON-05-020 TC-05-023 roomNo = 398 TCOV-05-002 TCON-05-002 charges = 2000.00		charges = 1999.99	TCOV-05-004	TCON-05-004
TCOV-05-014 TCON-05-014 TCOV-05-020 TCON-05-020 TC-05-023 roomNo = 398 TCOV-05-002 TCON-05-002 charges = 2000.00			TCOV-05-007	TCON-05-007
TCOV-05-020 TCON-05-020  TC-05-023 roomNo = 398 TCOV-05-002 TCON-05-002  charges = 2000.00			TCOV-05-009	TCON-05-009
TC-05-023 roomNo = 398 TCOV-05-002 TCON-05-002 charges = 2000.00			TCOV-05-014	TCON-05-014
charges = 2000.00			TCOV-05-020	TCON-05-020
charges = 2000.00	TC-05-023		TCOV-05-002	TCON-05-002
TCOV-05-004 TCON-05-004		charges = 2000.00	TCOV-05-004	TCON-05-004
TCOV-05-007 TCON-05-007			TCOV-05-007	TCON-05-007
TCOV-05-009 TCON-05-009			TCOV-05-009	TCON-05-009
TCOV-05-014 TCON-05-014			TCOV-05-014	TCON-05-014
TCOV-05-021 TCON-05-021			TCOV-05-021	TCON-05-021
TC-05-024 roomNo = 398 TCOV-05-002 TCON-05-002	TC-05-024		TCOV-05-002	TCON-05-002
charges = 2000.01 TCOV-05-004 TCON-05-004		charges = 2000.01	TCOV-05-004	TCON-05-004
TCOV-05-008 TCON-05-008			TCOV-05-008	TCON-05-008

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

		TCOV-05-010	TCON-05-010
		TCOV-05-014	TCON-05-014
		TCOV-05-022	TCON-05-022
TC-05-025	roomNo = 399	TCOV-05-002	TCON-05-002
	charges = 99.99	TCOV-05-004	TCON-05-004
		TCOV-05-006	TCON-05-006
		TCOV-05-010	TCON-05-010
		TCOV-05-015	TCON-05-015
		TCOV-05-017	TCON-05-017
TC-05-026	roomNo = 399	TCOV-05-002	TCON-05-002
	charges = 100.00	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-015	TCON-05-015
		TCOV-05-018	TCON-05-018
TC-05-027	roomNo = 399 charges = 100.01	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-015	TCON-05-015
		TCOV-05-019	TCON-05-019
TC-05-028	roomNo = 399	TCOV-05-002	TCON-05-002
	charges = 1999.99	TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID:



TCOV-05-009   TCON-05-009     TCOV-05-015   TCON-05-015     TCOV-05-020   TCON-05-020     TCOV-05-020   TCON-05-020     TCOV-05-020   TCON-05-002     TCOV-05-004   TCON-05-004     TCOV-05-007   TCON-05-007     TCOV-05-009   TCON-05-009     TCOV-05-015   TCON-05-009     TCOV-05-015   TCON-05-015     TCOV-05-021   TCON-05-021     TCOV-05-021   TCON-05-021     TCOV-05-021   TCON-05-021     TCOV-05-004   TCON-05-002     TCOV-05-004   TCON-05-002     TCOV-05-004   TCON-05-002     TCOV-05-004   TCON-05-004     TCOV-05-004   TCON-05-008     TCOV-05-008   TCON-05-008     TCOV-05-010   TCON-05-010     TCOV-05-015   TCON-05-015     TCOV-05-015   TCON-05-015     TCOV-05-016   TCON-05-005     TCOV-05-006   TCON-05-006     TCOV-05-010   TCON-05-010     TCOV-05-010   TCON-05-001     TCOV-05-010   TCON-05-001     TCOV-05-010   TCON-05-001     TCOV-05-001   TCON-05-001				
TC-05-029			TCOV-05-009	TCON-05-009
TC-05-029			TCOV-05-015	TCON-05-015
Charges = 2000.00			TCOV-05-020	TCON-05-020
TCOV-05-004 TCON-05-004  TCOV-05-007 TCON-05-007  TCOV-05-009 TCON-05-009  TCOV-05-015 TCON-05-015  TCOV-05-021 TCON-05-021  TCOV-05-021 TCON-05-021  TCOV-05-002 TCON-05-002  TCOV-05-004 TCON-05-002  TCOV-05-004 TCON-05-004  TCOV-05-008 TCON-05-008  TCOV-05-010 TCON-05-010  TCOV-05-015 TCON-05-015  TCOV-05-022 TCON-05-015  TCOV-05-022 TCON-05-022  TCOV-05-022 TCON-05-022  TCOV-05-022 TCON-05-022  TCOV-05-005 TCON-05-003  TCOV-05-006 TCON-05-006  TCOV-05-010 TCON-05-010	TC-05-029		TCOV-05-002	TCON-05-002
TCOV-05-009 TCON-05-009  TCOV-05-015 TCON-05-015  TCOV-05-021 TCON-05-021  TC-05-030 roomNo = 399 charges = 2000.01  TCOV-05-002 TCON-05-002  TCOV-05-004 TCON-05-004  TCOV-05-008 TCON-05-008  TCOV-05-010 TCON-05-010  TCOV-05-015 TCON-05-010  TCOV-05-022 TCON-05-015  TCOV-05-022 TCON-05-022  TCOV-05-003 TCON-05-003  TCOV-05-006 TCON-05-006  TCOV-05-010 TCON-05-010  TCOV-05-017 TCON-05-017  TC-05-032 roomNo = 400 charges = 100.00  TCOV-05-005 TCON-05-003  TCOV-05-005 TCON-05-005		charges = 2000.00	TCOV-05-004	TCON-05-004
TC-05-030			TCOV-05-007	TCON-05-007
TC-05-030			TCOV-05-009	TCON-05-009
TC-05-030			TCOV-05-015	TCON-05-015
Charges = 2000.01			TCOV-05-021	TCON-05-021
TCOV-05-004 TCON-05-004  TCOV-05-008 TCON-05-008  TCOV-05-010 TCON-05-010  TCOV-05-015 TCON-05-015  TCOV-05-022 TCON-05-022  TCOV-05-022 TCON-05-022  TCOV-05-031 TCOV-05-003 TCON-05-003  TCOV-05-006 TCON-05-006  TCOV-05-010 TCON-05-010  TCOV-05-010 TCON-05-010  TCOV-05-010 TCON-05-010  TCOV-05-017 TCON-05-017  TC-05-032 roomNo = 400 charges = 100.00  TCOV-05-003 TCON-05-003  TCOV-05-005 TCON-05-003	TC-05-030		TCOV-05-002	TCON-05-002
TCOV-05-010 TCON-05-010  TCOV-05-015 TCON-05-015  TCOV-05-022 TCON-05-022  TC-05-031 roomNo = 400		cnarges = 2000.01	TCOV-05-004	TCON-05-004
TCOV-05-015 TCON-05-015  TCOV-05-022 TCON-05-022  TC-05-031 roomNo = 400			TCOV-05-008	TCON-05-008
TC-05-031			TCOV-05-010	TCON-05-010
TC-05-031			TCOV-05-015	TCON-05-015
TCOV-05-005   TCON-05-005     TCOV-05-006   TCON-05-006     TCOV-05-010   TCON-05-010     TCOV-05-016   TCON-05-016     TCOV-05-017   TCON-05-017     TC-05-032   roomNo = 400			TCOV-05-022	TCON-05-022
TCOV-05-005 TCON-05-005  TCOV-05-006 TCON-05-006  TCOV-05-010 TCON-05-010  TCOV-05-016 TCON-05-016  TCOV-05-017 TCON-05-017  TC-05-032 roomNo = 400	TC-05-031		TCOV-05-003	TCON-05-003
TCOV-05-010 TCON-05-010 TCOV-05-016 TCON-05-016 TCOV-05-017 TCON-05-017 TC-05-032 roomNo = 400		charges = 99.99	TCOV-05-005	TCON-05-005
TCOV-05-016 TCON-05-016  TCOV-05-017 TCON-05-017  TC-05-032 roomNo = 400			TCOV-05-006	TCON-05-006
TCOV-05-017 TCON-05-017  TC-05-032 roomNo = 400			TCOV-05-010	TCON-05-010
TC-05-032 roomNo = 400 TCOV-05-003 TCON-05-003 TCON-05-005			TCOV-05-016	TCON-05-016
charges = 100.00 TCOV-05-005 TCON-05-005			TCOV-05-017	TCON-05-017
TCOV-05-005 TCON-05-005	TC-05-032		TCOV-05-003	TCON-05-003
TCOV-05-007 TCON-05-007		charges = 100.00	TCOV-05-005	TCON-05-005
			TCOV-05-007	TCON-05-007

Date:06/04/2024 Test Plan ID: KHS\_TDS\_1\_1.0.0

SINCE SON

		TCOV-05-009	TCON-05-009
		TCOV-05-016	TCON-05-016
		TCOV-05-018	TCON-05-018
TC-05-033	roomNo = 400	TCOV-05-003	TCON-05-003
	charges = 100.01	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-016	TCON-05-016
		TCOV-05-019	TCON-05-019
TC-05-034	roomNo = 400	TCOV-05-003	TCON-05-003
	charges = 1999.99	TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-016	TCON-05-016
		TCOV-05-020	TCON-05-020
TC-05-035	roomNo = 400 charges = 2000.00	TCOV-05-003	TCON-05-003
		TCOV-05-005	TCON-05-005
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-016	TCON-05-016
		TCOV-05-021	TCON-05-021
TC-05-036	roomNo = 400	TCOV-05-003	TCON-05-003
	charges = 2000.01	TCOV-05-005	TCON-05-005
		TCOV-05-008	TCON-05-008

Date:06/04/2024 Test Plan ID:

SINCE 202

		TCOV-05-010	TCON-05-010
		TCOV-05-016	TCON-05-016
		TCOV-05-022	TCON-05-022
TC-05-037	roomNo = 000 charges = "don't care"	TCOV-05-023	TCON-05-023
TC-05-038	roomNo = 400 charges = "don't care"	TCOV-05-025	TCON-05-025
TC-05-039	Add Room Main Flow	TCOV-05-002	TCON-05-002
		TCOV-05-004	TCON-05-004
		TCOV-05-007	TCON-05-007
		TCOV-05-009	TCON-05-009
		TCOV-05-028	TCON-05-028
TC-05-040	Add Room Alternate Flow -	TCOV-05-001	TCON-05-001
	Invalid roomNo	TCOV-05-003	TCON-05-003
		TCOV-05-005	TCON-05-005
		TCOV-05-029	TCON-05-029
TC-05-041	Add Room Alternate Flow -	TCOV-05-001	TCON-05-001
	Invalid charges	TCOV-05-004	TCON-05-004
		TCOV-05-006	TCON-05-006
		TCOV-05-008	TCON-05-008
		TCOV-05-010	TCON-05-010
		TCOV-05-030	TCON-05-030

Project Title: KHS Independent Verification and Validation		
Date:06/04/2024 Test Plan ID:		



### 2.4 Features Pass/Fail Criteria

KHS\_TDS\_1\_1.0.0

The system must satisfy the following criteria in order to pass:

## I. Add Room

- Pass Criteria: Successfully adds a new room to the system.
- Fail Criteria: Failure to add a room due to system error or incorrect data entry.

### II. Add Doctor

- Pass Criteria: Successfully adds a new doctor to the system.
- Fail Criteria: Failure to add a doctor due to system error or incomplete information.

## III. Add Patient

- Pass Criteria: Successfully adds a new patient to the system.
- Fail Criteria: Failure to add a patient due to system error or incomplete/incorrect information.

#### IV. Add Nurse

- Pass Criteria: Successfully adds a new nurse to the system.
- Fail Criteria: Failure to add a nurse due to system error or incomplete information.

### V. Add Medicine

- Pass Criteria: Successfully adds a new medicine to the system.
- Fail Criteria: Failure to add a medicine due to system error or incomplete information.

Project Title: KHS Independent Ver	rification and Validation
------------------------------------	---------------------------

Date:06/04/2024 Test Plan ID:

KHS\_TDS\_1\_1.0.0



## 2.5 Test Deliverables

The following documents will be generated by the system test group and will be delivered during the system testing.

## Test documentation:

- I. Test Plan
- II. Test Transmittal Report
- III. Test Design Specifications
- IV. Test Case Specifications
- V. Test Procedure Specifications
- VI. Test Logs
- VII. Test Incident Reports
- VIII. Test Summary Reports
- IX. Test Completion Report